

# Physics Lab Set-up

Thomas More College



# Physics Lab Set-up

Thomas More College

Joe Christensen  
Thomas More College

Credit for MathBookXML / PreTeXt format:

Robert A. Beezer

Latest update: February 26, 2018

**Edition:** Annual Edition 2017

**Website:** [TMC Physics](#)

© 2017–2018 J. Christensen

Permission is granted to copy, distribute and/or modify this document under the terms of the GNU Free Documentation License, Version 1.2 or any later version published by the Free Software Foundation; with no Invariant Sections, no Front-Cover Texts, and no Back-Cover Texts. A copy of the license is included in the appendix entitled “GNU Free Documentation License.” All trademarks™ are the registered® marks of their respective owners.

# Acknowledgements

I would like to acknowledge the following reviewers for their helpful comments and suggestions.

- Physics adjunct instructors who provided feedback about the labs
  - Tom Neal, Physics Adjunct
- Physics faculty who designed and contributed to the design of the lab experiments
  - Jack Wells
  - Dr. Wes Ryle
  - Dr. Jeremy Huber



# Download the PDF here

A large PDF of this document can be found at <http://physics.thomasmore.edu/labs/TMC-lab-setup.pdf>.



# Contents

|                                                                                    |            |
|------------------------------------------------------------------------------------|------------|
| <b>Acknowledgements</b>                                                            | <b>v</b>   |
| <b>Download the PDF here</b>                                                       | <b>vii</b> |
| <b>1 PHY 121L: Elements of Physics (algebra-based, fall)</b>                       | <b>1</b>   |
| 1.1 Meaningful Measurements . . . . .                                              | 1          |
| 1.2 Standard Deviation . . . . .                                                   | 2          |
| 1.3 Constant Acceleration . . . . .                                                | 3          |
| 1.4 Newton's 2 <sup>nd</sup> Law on a Linear Track with the Sonic Ranger . . . . . | 4          |
| 1.5 Dry Sliding Friction . . . . .                                                 | 5          |
| 1.6 Centripetal Force . . . . .                                                    | 9          |
| 1.7 Hooke's Law / Springs / Peer Review . . . . .                                  | 10         |
| 1.8 Pendulum / Peer Review . . . . .                                               | 13         |
| 1.9 Conservation of Energy . . . . .                                               | 17         |
| 1.10 Ballistic Pendulum . . . . .                                                  | 19         |
| 1.11 Human Forearm . . . . .                                                       | 22         |
| 1.12 Next Lab . . . . .                                                            | 27         |
| <b>2 PHY122L: Elements of Physics (algebra-based, spring)</b>                      | <b>29</b>  |
| 2.1 Skills Review . . . . .                                                        | 29         |
| 2.2 Archimedes' Principle . . . . .                                                | 29         |
| 2.3 Thermal Expansion . . . . .                                                    | 34         |
| 2.4 Specific and Latent Heat . . . . .                                             | 34         |
| 2.5 Properties of a Battery . . . . .                                              | 36         |
| 2.6 Simple Lens . . . . .                                                          | 38         |
| <b>3 NSC 220L: Anything Physics</b>                                                | <b>39</b>  |
| 3.1 Meaningful Measurements . . . . .                                              | 39         |
| 3.2 Measuring Motion . . . . .                                                     | 39         |
| 3.3 Constant Acceleration . . . . .                                                | 39         |
| 3.4 Newton's 2 <sup>nd</sup> Law on a Linear Track with the Sonic Ranger . . . . . | 40         |
| 3.5 Dry Sliding Friction . . . . .                                                 | 40         |
| 3.6 Hooke's Law . . . . .                                                          | 40         |
| 3.7 Pendulum / Peer Review . . . . .                                               | 40         |
| 3.8 Conservation of Energy . . . . .                                               | 40         |
| 3.9 Bouncing Springs . . . . .                                                     | 40         |
| <b>4 PHY141L: General Physics I (calculus-based, spring)</b>                       | <b>41</b>  |
| <b>5 PHY142L: General Physics II (calculus-based, spring)</b>                      | <b>43</b>  |
| 5.1 Ohm's Law (Ohmic material) . . . . .                                           | 43         |
| 5.2 Ohm's Law (non-Ohmic material) . . . . .                                       | 44         |
| 5.3 Mapping the Electric Field . . . . .                                           | 45         |
| 5.4 Resistor Bridge . . . . .                                                      | 47         |

|          |                                                              |           |
|----------|--------------------------------------------------------------|-----------|
| 5.5      | Magnetic Field Strength . . . . .                            | 48        |
| 5.6      | High Current Force . . . . .                                 | 48        |
| 5.7      | Magnet in a Tube . . . . .                                   | 50        |
| 5.8      | Capacitor and Variable Resistor . . . . .                    | 52        |
| <b>6</b> | <b>PHY241L: General Physics III (calculus-based, spring)</b> | <b>53</b> |
| 6.1      | Archimedes' Principle . . . . .                              | 53        |
| 6.2      | Standing Waves . . . . .                                     | 56        |
| 6.3      | Speed of Sound in a Tube . . . . .                           | 57        |
| 6.4      | Thermal Expansion . . . . .                                  | 59        |
| 6.5      | Boyle's Law . . . . .                                        | 60        |
| 6.6      | Specific Heat Capacity . . . . .                             | 62        |
| 6.7      | Newton's Law of Cooling . . . . .                            | 64        |
| 6.8      | Snell's Law of Refraction . . . . .                          | 65        |
| 6.9      | Reflection and Refraction . . . . .                          | 68        |
| 6.10     | Radiation . . . . .                                          | 70        |
| 6.11     | Interferometer . . . . .                                     | 70        |


# Class 1

## PHY 121L: Elements of Physics (algebra-based, fall)


### 1.1 Meaningful Measurements

| Location                                         | Equipment                                          | Notes                                                                                                                                                                               |
|--------------------------------------------------|----------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b><i>For Each Lab Station</i></b>               |                                                    |                                                                                                                                                                                     |
| AA13                                             | 1 Vernier caliper ( <a href="#">Figure 1.1.2</a> ) | Ask faculty if they want digit or analog before and after lab, verify digital calipers are turned off                                                                               |
| AA14                                             | 1 micrometer ( <a href="#">Figure 1.1.2</a> )      | before and after lab verify jaws are not tight                                                                                                                                      |
| AF34-14                                          | 1 metric ruler                                     | could be a 1-foot ruler or a 0.25-meter stick                                                                                                                                       |
| AF55-19                                          | 1 tall, skinny, graduated cylinder                 | at least 2 of the 3 objects should fit inside the graduated cylinder                                                                                                                |
| AE21                                             | string                                             | sufficient string to tie onto the objects in order to immerse them in the cylinder and retrieve them                                                                                |
| AE82                                             | 3 objects to measure                               | Ask faculty which objects they want (smooth sphere, rough sphere, cube, block of metal, irregular shape, etc) <b><i>These objects should fit inside the graduated cylinder.</i></b> |
| <b><i>At the front for students to share</i></b> |                                                    |                                                                                                                                                                                     |
| S224                                             | at least 1 digital scale                           | the available scale(s) should be able to weigh the chosen objects                                                                                                                   |
| AF36-4                                           |                                                    |                                                                                                                                                                                     |

**Table 1.1.1:** Equipment Needed: Meaningful Measurements

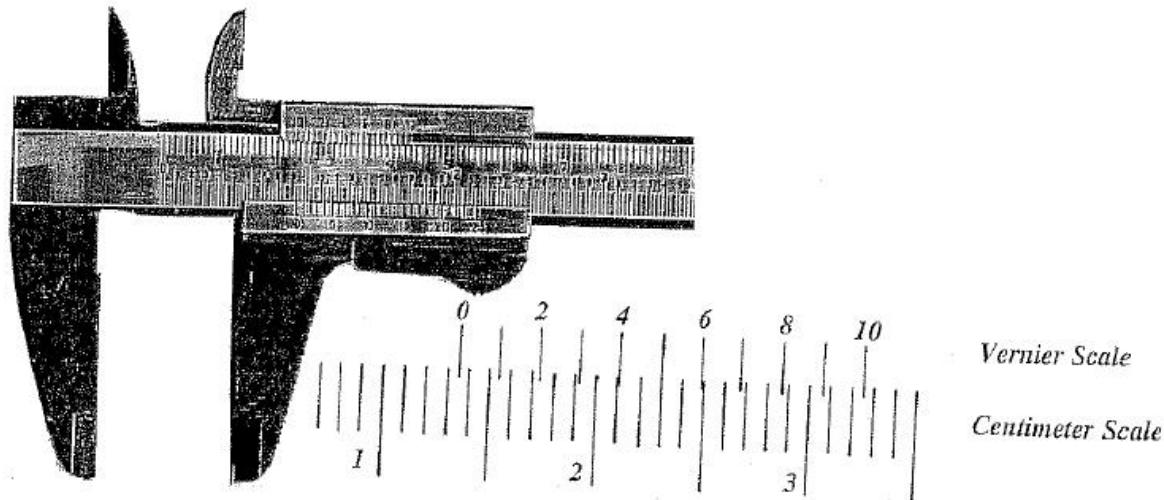

**Figure 1.1.2:** The Vernier Caliper

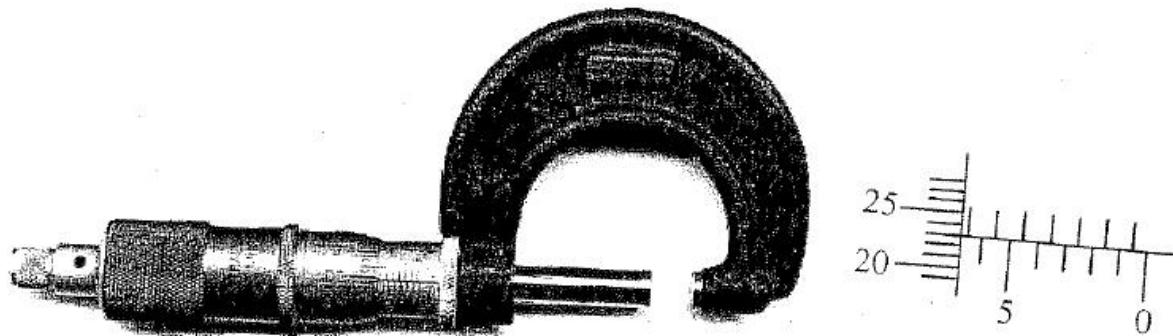

**Figure 1.1.3:** The Micrometer Caliper

(Updated: October 12, 2017)

A digital version of the lab should be found at <http://physics.thomasmore.edu/Labs/121/c-meaningful-measurements.html>

A PDF version of the write-up might be found at [Measurement.pdf \(291 kB\)](#)

## 1.2 Standard Deviation

| Location                                  | Equipment                     | Notes                                                          |
|-------------------------------------------|-------------------------------|----------------------------------------------------------------|
| <i>At the front for students to share</i> |                               |                                                                |
| AL32-5                                    | several boxes/bags of pennies | in sufficient number for each student to have up to 50 pennies |

**Table 1.2.1:** Equipment Needed: Standard Deviation

(Updated: September 6, 2017)

A digital version of the lab should be found at <http://physics.thomasmore.edu/Labs/121/c-standard-deviation.html>

A PDF version of the write-up might be found at [StDev.pdf \(232 kB\)](#)

## 1.3 Constant Acceleration

| Location                                  | Equipment                                        | Notes                                                                                                                                                                                                                                                                   |
|-------------------------------------------|--------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>For Each Lab Station</b>               |                                                  |                                                                                                                                                                                                                                                                         |
| S224                                      | 1 motion sensor (same as “sonic ranger”)         | Should have a black-yellow plug                                                                                                                                                                                                                                         |
| AF12 (tube)                               | 1 track                                          | .                                                                                                                                                                                                                                                                       |
| AF15                                      | Wood Squares                                     | probably 2-3, used to prop up one end of track                                                                                                                                                                                                                          |
| AF22-2-13                                 | 1 cart with “sail”                               | these are in a large box labelled “DYNAMIC CARTS AF22-2-13”                                                                                                                                                                                                             |
| AF34-14                                   | ruler                                            | used to level the track                                                                                                                                                                                                                                                 |
| either                                    |                                                  |                                                                                                                                                                                                                                                                         |
| AF35<br>(shelf) or<br>AE82<br>(drawer)    | 1 metal ball (any size)                          | used to level the track                                                                                                                                                                                                                                                 |
| -                                         | Pasco                                            | Computer                                                                                                                                                                                                                                                                |
| <b>At the front for students to share</b> |                                                  |                                                                                                                                                                                                                                                                         |
| AA41 or<br>AA42<br>(drawers)              | 1 gravity protractor                             | This is the large yellow protractor                                                                                                                                                                                                                                     |
| AF44                                      | Pendulum bob to hang with the gravity-protractor | Used to level the track. (There are enough protractors to give these to everybody, but it is not clear if there are enough pendula strings and bobs for everybody.)<br><i>This really only needs to be a short string with a washer or other small mass tied to it.</i> |

**Table 1.3.1:** Equipment Needed: [Constant Acceleration](#)

(Updated: October 12, 2017)


A digital version of the lab should be found at <http://physics.thomasmore.edu/Labs/121/c-acceleration.html>

A PDF version might be found at [Acceleration.pdf](#)

## 1.4 Newton's 2<sup>nd</sup> Law on a Linear Track with the Sonic Ranger

| Location                                         | Equipment                                        | Notes                                                                                                                                                                                                                                                                    |
|--------------------------------------------------|--------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>For Each Lab Station</b>                      |                                                  |                                                                                                                                                                                                                                                                          |
| S224                                             | 1 motion sensor (AKA "sonic ranger")             | Should have a black-yellow plug                                                                                                                                                                                                                                          |
| AF12 (tube)                                      | 1 track                                          | .                                                                                                                                                                                                                                                                        |
| AF15                                             | 2-3 small wooden squares                         | used to level the track                                                                                                                                                                                                                                                  |
| AF22-2-13                                        | 1 cart                                           | these are in a large box labelled "DYNAMIC CARTS AF22-2-13"                                                                                                                                                                                                              |
| AF22-2-13                                        | 1 wooden cart-block                              | these are in a large box labelled "DYNAMIC CARTS AF22-2-13"                                                                                                                                                                                                              |
| AF22-2-13                                        | light plastic bucket                             | These might already be attached to the string in a large box labelled "DYNAMIC CARTS AF22-2-13"                                                                                                                                                                          |
| (attached<br>to plastic<br>bucket?)              | string                                           | There should be pre-cut string that is long enough to reach from the cart, over the pulley and to a hanging mass. About one meter long. This is likely in a large box labelled "DYNAMIC CARTS AF22-2-13"                                                                 |
| AF22-2-13                                        | 1 pulley                                         | these are in a large box labelled "DYNAMIC CARTS AF22-2-13"                                                                                                                                                                                                              |
| AF34-14                                          | ruler                                            | used to level the string                                                                                                                                                                                                                                                 |
| either<br>AF35<br>(shelf) or<br>AE82<br>(drawer) | 1 metal ball (any size)                          | used to check the level of the track                                                                                                                                                                                                                                     |
| AF44                                             | larger weights                                   | These are to ride the cart. <b>Check with instructor: EITHER</b> an assortment of 100-500 gram, cylindrical masses <b>OR</b> 2 black rectangular masses that fit into the cart (like the 1-wooden block above).                                                          |
| AF44 or<br>AL32-5                                | tiny weights                                     | These are to transfer between the cart and the basket. <b>Check with instructor: EITHER</b> an assortment of 7-10 very small masses (2-5 grams) <b>OR</b> 10 pennies.                                                                                                    |
| -                                                | Pasco                                            | Computer                                                                                                                                                                                                                                                                 |
| <b>At the front for students to share</b>        |                                                  |                                                                                                                                                                                                                                                                          |
| S224                                             | functioning digital scales                       | (please verify that these function and are set to metric)                                                                                                                                                                                                                |
| AA41,<br>AA42                                    | Gravity Protractor                               | Used to level the track. (There are enough to give these to everybody, but it is not clear if there are enough pendulua strings and bobs for everybody.)                                                                                                                 |
| AF44                                             | Pendulum bob to hang with the gravity-protractor | Used to level the track. (There are enough protractors to give these to everybody, but it is not clear if there are enough pendulua strings and bobs for everybody.)<br><i>This really only needs to be a short string with a washer or other small mass tied to it.</i> |

**Table 1.4.1:** Equipment Needed: Newton's 2<sup>nd</sup> Law on a Linear Track with the Sonic Ranger

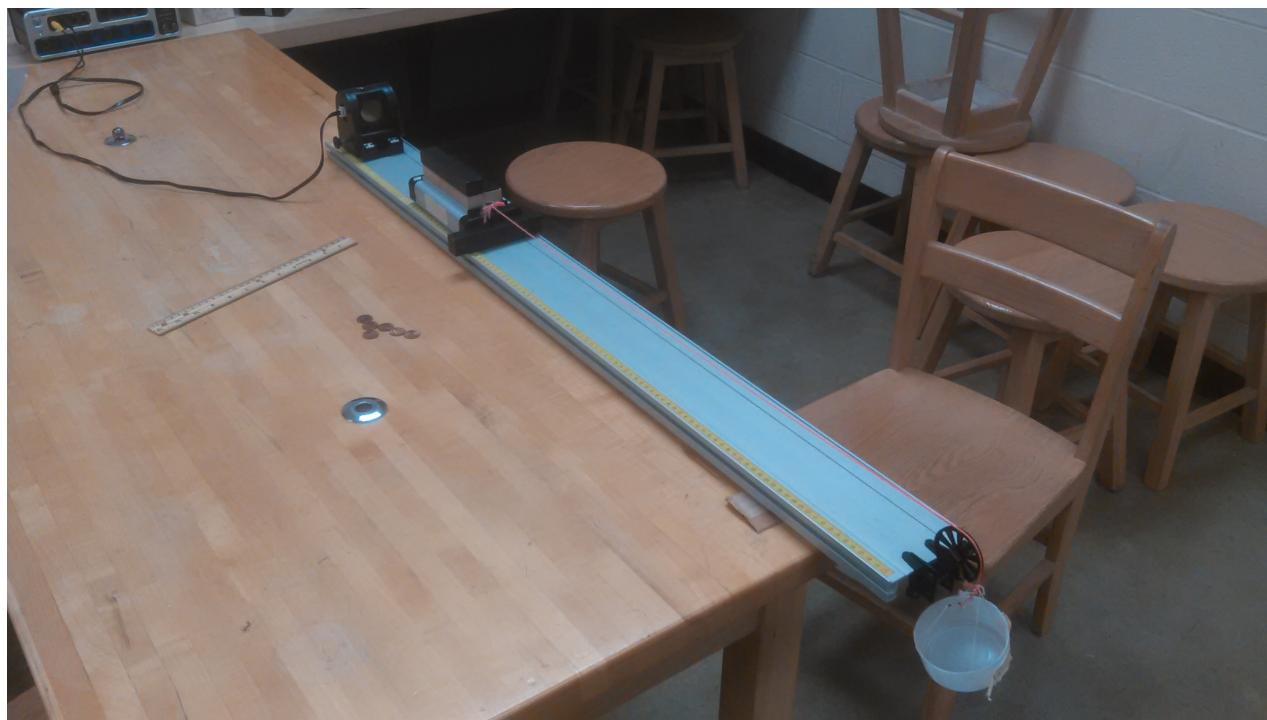

**Figure 1.4.2:** The equipment setup.

(Updated: October 12, 2017)

A digital version of the lab should be found at <http://physics.thomasmore.edu/Labs/121/c-Newton.html>

A PDF version might be found at [Newton.pdf](#)

## 1.5 Dry Sliding Friction


Most of this equipment is the same as [Newton's 2<sup>nd</sup> Law on a Linear Track with the Sonic Ranger](#), but it is not exactly the same.

| Location                                         | Equipment                                        | Notes                                                                                                                                                                                                                                                                        |
|--------------------------------------------------|--------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>For Each Lab Station</b>                      |                                                  |                                                                                                                                                                                                                                                                              |
| S224                                             | 1 motion sensor (AKA “sonic ranger”)             | Should have a black-yellow plug                                                                                                                                                                                                                                              |
| S224                                             | 1 force transducer                               | Should be the thing shaped to sit on top of the wooden cart-block and it has a loop-screw at one side to tie a string to                                                                                                                                                     |
| AF12 (tube)                                      | 1 track                                          | .                                                                                                                                                                                                                                                                            |
| AF22-2-13                                        | 1 cart                                           | these are in a large box labelled “DYNAMIC CARTS AF22-2-13”                                                                                                                                                                                                                  |
| AF22-2-13                                        | 1 wooden cart-block                              | these are in a large box labelled “DYNAMIC CARTS AF22-2-13”                                                                                                                                                                                                                  |
| AF22-2-13                                        | light plastic bucket                             | These might already be attached to the string in a large box labelled “DYNAMIC CARTS AF22-2-13”                                                                                                                                                                              |
| (attached to plastic bucket?)                    | string                                           | There should be pre-cut string that is long enough to reach from the cart, over the pulley and to a hanging mass. About one meter long. This is likely in a large box labelled “DYNAMIC CARTS AF22-2-13”                                                                     |
| AF22-2-13                                        | 1 pulley                                         | these are in a large box labelled “DYNAMIC CARTS AF22-2-13”                                                                                                                                                                                                                  |
| AF34-14                                          | ruler                                            | used to level the string                                                                                                                                                                                                                                                     |
| either<br>AF35<br>(shelf) or<br>AE82<br>(drawer) | 1 metal ball (any size)                          | used to check the level of the track                                                                                                                                                                                                                                         |
| AF15                                             | 2-3 small wooden squares                         | used to level the track                                                                                                                                                                                                                                                      |
| AF44                                             | larger weights                                   | These are to ride the cart. <b>Check with instructor: EITHER</b> cylindrical masses (2 500-g, 4 200-g, 2 100-g) <b>OR</b> 2 black rectangular masses that fit into the cart (like the 1-wooden block above).                                                                 |
| AF44 or<br>AL32-5                                | tiny weights                                     | These are to transfer between the cart and the basket. <b>Check with instructor: EITHER</b> an assortment of 7-10 very small masses (2-5 grams) <b>OR</b> 10 pennies.                                                                                                        |
| AE21                                             | string                                           | This string is for pulling the force transducer, not the same as the string that is attached to the plastic bucket                                                                                                                                                           |
| AG15                                             | 1 thermal glove                                  | <b>Check with the instructor.</b> This is only used as a soft glove to catch the cart in case it accelerates too quickly. Not all instructors want this.                                                                                                                     |
| -                                                | Pasco                                            | Computer                                                                                                                                                                                                                                                                     |
| <b>At the front for students to share</b>        |                                                  |                                                                                                                                                                                                                                                                              |
| S224                                             | functioning digital scales                       | (please verify that these function and are set to metric)                                                                                                                                                                                                                    |
| AA41,<br>AA42                                    | Gravity Protractor                               | Used to level the track. (There are enough to give these to everybody, but it is not clear if there are enough pendulua strings and bobs for everybody.)                                                                                                                     |
| AF44                                             | Pendulum bob to hang with the gravity-protractor | Used to level the track. (There are enough protractors to give these to everybody, but it is not clear if there are enough pendulua strings and bobs for everybody.)<br><br><i>This really only needs to be a short string with a washer or other small mass tied to it.</i> |

Table 1.5.1: Equipment Needed: Dry Sliding Friction



**Figure 1.5.2:** The equipment set up.

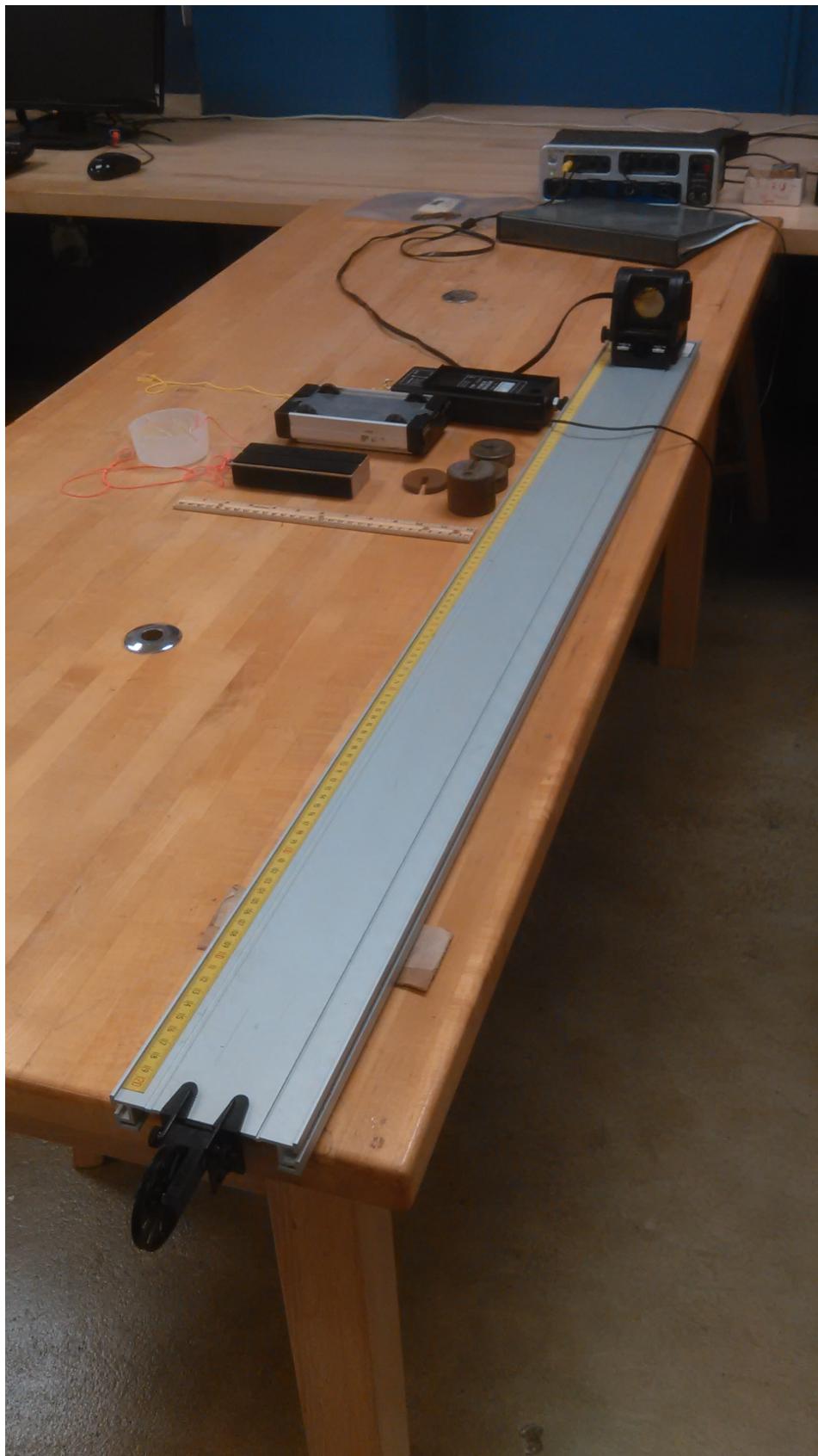

**Figure 1.5.3:** The equipment setup.

(Updated: October 12, 2017)

A digital version of the lab should be found at <http://physics.thomasmore.edu/Labs/121/c-friction.html>

A PDF version might be found at [friction.pdf](#)

## 1.6 Centripetal Force

This describes the old equipment. New equipment (that plugs into the PASCO interface) has been purchased, but that equipment list has not been itemized yet.

| Location                    | Equipment                                               | Notes                                                                                                                                                   |
|-----------------------------|---------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>For Each Lab Station</b> |                                                         |                                                                                                                                                         |
| AA13                        | 1 caliper                                               | <b>Ask instructor</b> if they want digital or analog (they might not care)                                                                              |
| AA12                        | 1 stop watch                                            | Verify that they work and have good batteries.<br><b>You might need to purchase new batteries for department!!</b>                                      |
| AD71                        | 1 C-clamp                                               | <b>Important:</b> There are small pieces of cardboard available to place between the clamp and the lab table so that you do not add dents to the table! |
| end of AD                   | 1 thick metal rod                                       | Mount horizontally off the table to hang the (spinning) carriage                                                                                        |
| AF14-1-6                    | 1 centripetal force apparatus                           | (with card?)                                                                                                                                            |
| AF46                        | 1 0.5-kg hanger                                         | .                                                                                                                                                       |
| AF46                        | Large cylindrical masses: 1 2-kg,<br>2 1-kg, 1 0.5-kg   | .                                                                                                                                                       |
| AF44                        | Additional small cylindrical<br>masses of various sizes | include at least one of each: 200-g, 100-g, 50-g, 20-g                                                                                                  |

**Table 1.6.1:** Equipment Needed: Centripetal Force



Figure 1.6.2: The equipment setup.

(Updated: October 12, 2017)

A digital version of the lab should be found at <http://physics.thomasmore.edu/Labs/121/c-centripetal.html>

A PDF version might be found at [centripetal.pdf](#)

## 1.7 Hooke's Law / Springs / Peer Review

**Warning 1.7.1.** This lab might be set-up by itself or it might be paired with [Pendulum / Peer Review](#). **Check with your instructor** to see if this is part of the “Peer Review” setup. If so, then discuss the setup. The instructor might want both Springs and Pendulum in both the morning and afternoon labs (in which case you need to verify how many setups of each lab and which tables), or the instructor might want one lab setup in the morning and a different lab setup in the afternoon.

| Location                                  | Equipment                    | Notes                                                                                                                                                                                                     |
|-------------------------------------------|------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>For Each Lab Station</b>               |                              |                                                                                                                                                                                                           |
| AD41                                      | 1 2-slot-perpendicular clamp | This is used to connect the vertical thick rod to the horizontal small rod                                                                                                                                |
| AD71                                      | 1 C-clamp                    | <b>Important:</b> There are small pieces of cardboard available to place between the clamp and the lab table so that you do not add dents to the table!                                                   |
| AD82                                      | 1 clamp                      | This is used to hold the track in place. Note the size of rod that fits through the hole on this clamp. <b>CAUTION:</b> Do not over-tighten this clamp. <i>The track should not be bent by the clamp!</i> |
| end of AD                                 | 1 thick metal rod            | Mount vertically off the table                                                                                                                                                                            |
| end of AD                                 | 1 small rod                  | This is used to connect the clamp to the thick rod                                                                                                                                                        |
| AF12                                      | 1 track                      | This is used to measure the elongation of the spring.                                                                                                                                                     |
| AF44                                      | hangers and various masses   | each table should have a similar distribution of large and small mass sizes ( <b>Check with instructor:</b> Your instructor might prefer to use the AF34 box of black hanging cylinders.)                 |
| <b>At the front for students to share</b> |                              |                                                                                                                                                                                                           |
| AF43                                      | 2 boxes of colored springs   | Students will be allowed to select their spring. <b>Check with instructor:</b> Your instructor might want specific springs distributed to the tables                                                      |

**Table 1.7.2:** Equipment Needed: [Hooke's Law / Springs / Peer Review](#)



**Figure 1.7.3:** The equipment setup version “A”. Your instructor might want the track setup in a different way.



**Figure 1.7.4:** The equipment setup version “A”. Your instructor might want the track setup in a different way.

(Updated: October 12, 2017)

A digital version of the lab should be found at <http://physics.thomasmore.edu/Labs/121/c-springs.html>

A PDF version might be found at [springs.pdf](#)

## 1.8 Pendulum / Peer Review


**Warning 1.8.1.** This lab might be set-up by itself or it might be paired with [Hooke’s Law / Springs / Peer Review](#). **Check with your instructor** to see if this is part of the “Peer Review” setup. If so, then discuss the setup. The instructor might want both Springs and Pendulum in both the morning and afternoon labs (in which case you need to verify how many setups of each lab and which tables), or the instructor might want one lab setup in the morning and a different lab setup in the afternoon.

| Location                    | Equipment                     | Notes                                                                                                                                                                                                                                               |
|-----------------------------|-------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>For Each Lab Station</b> |                               |                                                                                                                                                                                                                                                     |
| AA??,<br>maybe<br>AA41?     | 1 plastic protractor          | To measure the amplitude of the pendulum swing.                                                                                                                                                                                                     |
| AA12                        | 1 stop watch                  | Verify that they work and have good batteries.<br><br><i>You might need to purchase new batteries for department!!</i>                                                                                                                              |
| AD21                        | 1 3-notch black “arm” support | See picture below for how to mount the string inside the clamp.<br>Do not just wrap it around the screw; clamp it inside the metal piece, which you should note has a nub holding it in place. (Take the metal square off, you’ll see what I mean.) |
| AD71                        | 1 C-clamp                     | <b>Important:</b> There are small pieces of cardboard available to place between the clamp and the lab table so that you do not add dents to the table!                                                                                             |
| end of AD                   | 1 thick metal rod             | Mount vertically off the table                                                                                                                                                                                                                      |
| AF44                        | hangers and various masses    | each table should have a similar distribution of large and small mass sizes ( <b>Check with instructor:</b> Your instructor might prefer to use the AF34 box of black hanging cylinders.)                                                           |

**Table 1.8.2:** Equipment Needed: [Pendulum / Peer Review](#)



**Figure 1.8.3:** The equipment setup.

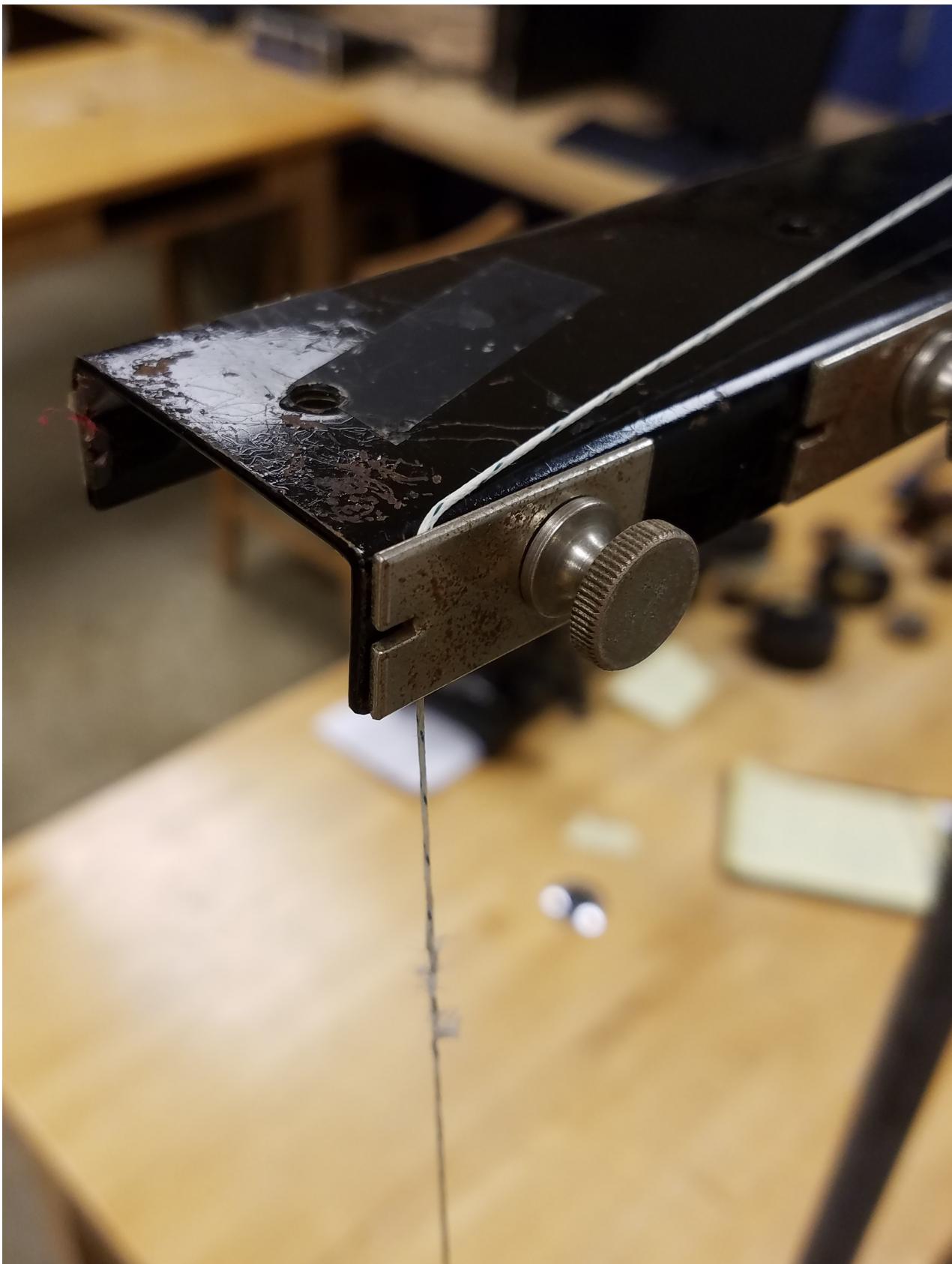

**Figure 1.8.4:** The equipment setup.

(Updated: October 12, 2017)

A digital version of the lab should be found at <http://physics.thomasmore.edu/Labs/121/c-pendulum.html>

A PDF version might be found at [pendulum.pdf](#)

## 1.9 Conservation of Energy

**Check with the instructor:** There are two versions of this lab. One is designed for a single week; the other is designed to do over two-weeks. All of the equipment used in the first week would also be used in the second week, but if you are doing the 2-week version, there are two additional items needed during the second week.

**Check with the instructor:** This lab is written with the intention of having students decide how to measure certain variables. This arrangement of equipment assumes that some equipment will be at the front of the room available for students to select; however, your instructor might prefer to have all of it at each lab station. Please confirm with your specific instructor.

| Location                                                 | Equipment                                        | Notes                                                                                                                                                                                                                                                                        |
|----------------------------------------------------------|--------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>For Each Lab Station</b>                              |                                                  |                                                                                                                                                                                                                                                                              |
| AF12 (tube)                                              | 1 track                                          | .                                                                                                                                                                                                                                                                            |
| AF22-2-13                                                | 1 cart                                           | these are in a large box labelled "DYNAMIC CARTS AF22-2-13"                                                                                                                                                                                                                  |
| AF22-2-13                                                | light plastic bucket                             | These might already be attached to the string in a large box labelled "DYNAMIC CARTS AF22-2-13"                                                                                                                                                                              |
| (attached<br>to plastic<br>bucket?)                      | string                                           | There should be pre-cut string that is long enough to reach from the cart, over the pulley and to a hanging mass. About one meter long. This is likely in a large box labelled "DYNAMIC CARTS AF22-2-13"                                                                     |
| AF22-2-13                                                | 1 pulley                                         | these are in a large box labelled "DYNAMIC CARTS AF22-2-13"                                                                                                                                                                                                                  |
| AF22<br>either<br>AF35<br>(shelf) or<br>AE82<br>(drawer) | C-shaped magnetic for the track                  | used in the second part of the lab (it will replace the pulley)                                                                                                                                                                                                              |
| AF15                                                     | 2-3 small wooden squares                         | used to level the track                                                                                                                                                                                                                                                      |
| -                                                        | Pasco                                            | Computer                                                                                                                                                                                                                                                                     |
| <b>At the front for students to share</b>                |                                                  |                                                                                                                                                                                                                                                                              |
| S224                                                     | functioning digital scales                       | (please verify that these function and are set to metric)                                                                                                                                                                                                                    |
| S224                                                     | 1 motion sensor (AKA "sonic ranger")             | Should have a black-yellow plug                                                                                                                                                                                                                                              |
| S224                                                     | 1 force transducer                               | <b>This is only needed in the second week</b> (if you are doing a second week)                                                                                                                                                                                               |
| AD91                                                     | meterstick                                       | used to measure the height of the hanging basket                                                                                                                                                                                                                             |
| AF34-14                                                  | ruler                                            | used to level the string                                                                                                                                                                                                                                                     |
| AA41,<br>AA42                                            | Gravity Protractor                               | Used to level the track. (There are enough to give these to everybody, but it is not clear if there are enough pendulua strings and bobs for everybody.)                                                                                                                     |
| AF44                                                     | Pendulum bob to hang with the gravity-protractor | Used to level the track. (There are enough protractors to give these to everybody, but it is not clear if there are enough pendulua strings and bobs for everybody.)<br><br><i>This really only needs to be a short string with a washer or other small mass tied to it.</i> |
| AF44                                                     | larger weights                                   | These are to ride the cart. A few cylindrical masses in the 200-g to 500-g range.                                                                                                                                                                                            |
| AF44                                                     | larger weights                                   | These are to ride the cart. Two black rectangular masses that fit into the cart.                                                                                                                                                                                             |
| AF44 or<br>AL32-5                                        | tiny weights                                     | They will only need one or two to accelerate the basket. Provide an assortment of 7-10 very small masses (2-5 grams). (Pennies would work, too.)                                                                                                                             |

**Table 1.9.1:** Equipment Needed: [Conservation of Energy](#)



**Figure 1.9.2:** The equipment setup.

(Updated: October 25, 2017)

A digital version of the 1-week version of the lab should be found at <http://physics.thomasmore.edu/Labs/121/c-energy-1.html>

A PDF version might be found at [energy-1.pdf](#)

A digital version of the 2-week version of the lab should be found at <http://physics.thomasmore.edu/Labs/121/c-energy-2.html>

A PDF version might be found at [energy-2.pdf](#)

## 1.10 Ballistic Pendulum

There is **additional work** for you to do beyond simply setting out the equipment this week.

- You will need to adjust the screws at the top of the pendulum arm.
  - While making this adjustment, you should remove the pendulum arm and verify the mass that is written on the masking tape. If that number is wrong, please update the tape and date it.
  - When you are making the adjustments, you will need to shoot the gun multiple times. Please note the following tips for cocking and shooting the gun.
    - Put the ball on the gun tip to cushion your hand when cocking the gun.
    - Brace your “other” hand on the back of the trigger-handle while cocking the gun.
    - The gun, when shot, has enough kick to bump the whole apparatus back a little, which can cause the pendulum arm to swing before the ball is caught, which might cause the ball to not be caught. Clamping the apparatus to the table will minimize this.
    - If you clamp the apparatus to the table, the clamp needs to be out of the path of the gun – this affects the way the tables are arranged in the room.
    - Be aware of what the gun might hit if it misses the pendulum bob. Do not shoot it towards people or expensive equipment!
  - Your goal while adjusting these screws is to simultaneously allow

1. the gun to shoot *straight* into the pendulum bob AND
  2. the pendulum to swing back and catch on the notches of the angled platform.
- You *might* need to trade the arms between the apparatuses. The steel machines are interchangeable and the red machines are interchangeable, but you cannot interchange the arms between red and steel.
- You will need to rearrange the lab tables so that students can shoot the gun giving the ball a 2-dimensional trajectory off the end and onto the floor. Please consider the following: (Assuming the lab room door is on the ***north*** side of the room)
    - The northside tables can shoot northwards. Be sure that the middle table is not aimed at the glass doors. You might have to move it slightly to the east (front of the room). The apparatuses should be clamped to the north east corner of the table. They do not need to be ALL the way north on the table, but if they are too far back from the edge, they might hit the table.
    - The three south tables will need to be moved. You have some options. In each case, you have to remember to leave room for the lab students to move around without crossing paths and access both their gun and the box that catches the ball being shot.
      - The south east and south west tables might be aligned long-ways east-west (rather than long-ways north-south as they are now). In this case, the south west table (apparatus clamped on the NW corner) shoots west (back of the room) and the south east table (apparatus clamped on the SE corner) shoots east (towards the white board). These should be pulled to the north-south middle of the room so that the south middle table can be aligned long-ways east-west to the south of these tables and can shoot either east (apparatus clamped on the SE corner) or west (apparatus clamped on the NW corner).
      - The south east and south west tables might be aligned long-ways north-south (as they are now); but they should be squeezed to abut the south middle table. In this case, the tables should be slightly off-set in order to allow students to access their guns and the location their gun fires to, without interfering with their neighbor. In this case, the south west table (apparatus clamped on the NW corner) shoots west (back of the room) and the south east table (apparatus clamped on the SE corner) shoots east (towards the white board). The south-west and south-east tables should be moved slightly towards the north-south middle of the room so that the south middle table can stick out a bit to the south of these tables. The south middle table should shoot east (towards the white board) specifically so that it can be clamped to the south-east corner of the table.
      - It is possible to arrange the room so that the southern table shoot to the south, but then the students have to crawl under the counter and they are shooting towards the computers along the south wall. I recommend against this option.
      - When you move the tables, ***be aware*** that the tables have electrical outlets in them and ***all have plugs*** that are in the floor or in the wall. ***Be careful not to pull the plugs out of the floor*** when moving the tables. If you get a small screw driver, then you can unwind some extra cable and relocate the stopper that keeps the cable from being re-wound. (See [Figure 3](#) below.)
    - All tables will need to be adjusted so that there is sufficient room to shoot and land far enough from the wall that a box can be placed between the landing spot and the wall. If the landing spot is too close to the box, then the ball *could* bounce off the wall and back onto the carbon paper, which makes the lab more difficult than necessary. Please leave enough room between the landing spot and the box that protects the wall that the ball does not immediately drop back onto the carbon paper that will be positioned at the landing site.

| Location                                  | Equipment                                           | Notes                                                                                                                                                                                                                                                                                                                                                                                          |
|-------------------------------------------|-----------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>For Each Lab Station</b>               |                                                     |                                                                                                                                                                                                                                                                                                                                                                                                |
| AF26                                      | Ballistic Pendulum apparatus                        | One is missing a foot and uses a black rubber cork as a stabilizer. Several have difficulty with keeping the pendulum portion fixed in place. You are expected to check the maintenance of this equipment early enough before lab that you (or possibly the instructor) can fix any flaws in the equipment before the lab. <b><i>Do not wait until the last minute to set this lab up!</i></b> |
| AD91,<br>AF34                             | 1 1/2-meter stick, 1 1-meter stick, 1 2-meter stick | Note they might want more metersticks.                                                                                                                                                                                                                                                                                                                                                         |
| AD72<br><b>AD43</b>                       | 1 C-clamp                                           | The AD72 clamps will work for the red ballistic pendulums, because they are short. The steel ones, though, are too tall and need the AD43 C-clamps. (I might be wrong about the AD43 bin number; please correct me so I can update this file.)                                                                                                                                                 |
|                                           |                                                     | <b>Important:</b> There are small pieces of cardboard available to place between the clamp and the lab table so that you do not add dents to the table!                                                                                                                                                                                                                                        |
| ?                                         | pendulum bob                                        | When you find this, tell Dr. Christensen the bin-number so he can update this page.                                                                                                                                                                                                                                                                                                            |
| ?                                         | carbon paper                                        | When you find this, tell Dr. Christensen the bin-number so he can update this page.                                                                                                                                                                                                                                                                                                            |
| take from<br>the printer<br>in S224A      | 5 sheets of printer paper                           |                                                                                                                                                                                                                                                                                                                                                                                                |
| <b>At the front for students to share</b> |                                                     |                                                                                                                                                                                                                                                                                                                                                                                                |
| ?                                         | bubble level                                        | When you find this, tell Dr. Christensen the bin-number so he can update this page.                                                                                                                                                                                                                                                                                                            |
| ?                                         | 2 rolls of masking tape                             | When you find this, tell Dr. Christensen the bin-number so he can update this page.                                                                                                                                                                                                                                                                                                            |

**Table 1.10.1:** Equipment Needed: Ballistic Pendulum



Figure 1.10.2: The equipment setup.



Figure 1.10.3: When rearranging the tables, but aware that under each table, there is a cord plugged into the floor or the wall. These can be unwound and, with a screwdriver, the stopper can be moved so that a longer stretch of cord is extended.

(Updated: November 15, 2017)

A digital version of the lab should be found at <http://physics.thomasmore.edu/Labs/121/c-ballistic.html>  
A PDF version might be found at [ballistic.pdf](#)

## 1.11 Human Forearm

For the instructor: We do **not** have permission to reproduce this image, but it might be useful for you to have in your head when explaining the lab.

| Location                                  | Equipment                     | Notes                                                                                                                                                                                                                                                              |
|-------------------------------------------|-------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>For Each Lab Station</b>               |                               |                                                                                                                                                                                                                                                                    |
| AD21                                      | 3 clamps for the meter-sticks | These are the little aluminum things that slide onto a meter-stick and have a metal piece that a mass-hanger can hang off of. If the instructor wants you (instead of the students to put these on, notice that two hang down and one hangs up – see image below.) |
| AD41                                      | 1 2-slot-perpendicular clamp  | This is used to connect the vertical thick rod to the horizontal small rod                                                                                                                                                                                         |
| AD71                                      | 1 C-clamp                     | <b>Important:</b> There are small pieces of cardboard available to place between the clamp and the lab table so that you do not add dents to the table!                                                                                                            |
| AD82                                      | 1 clamp                       | This is used to hold the track in place. Note the size of rod that fits through the hole on this clamp. <b>CAUTION:</b> Do not over-tighten this clamp. <i>The track should not be bent by the clamp!</i>                                                          |
| AD91                                      | 1 1-meter stick               | .                                                                                                                                                                                                                                                                  |
| end of AD                                 | 1 thick metal rod             | Mount vertically off the table                                                                                                                                                                                                                                     |
| end of AD                                 | 1 small rod                   | This is used to connect the clamp to the thick rod                                                                                                                                                                                                                 |
| AF12                                      | 1 track                       | This is used to measure the elongation of the spring.                                                                                                                                                                                                              |
| AF22-2-13                                 | light plastic bucket          | These likely have string attached and are in a large box labelled "DYNAMIC CARTS AF22-2-13"                                                                                                                                                                        |
| AF34                                      | 1 1/2-meter stick             | (might need a 1/2-meter stick)                                                                                                                                                                                                                                     |
| AF43                                      | 1 spring                      | (any color)                                                                                                                                                                                                                                                        |
| AF43                                      | 1 small rectangular mirror    | These are for minimizing parallax.                                                                                                                                                                                                                                 |
| AF44                                      | hangers and various masses    | The hangers might be in AF43, but their home is AF44.                                                                                                                                                                                                              |
| <b>At the front for instructor to use</b> |                               |                                                                                                                                                                                                                                                                    |
| AK32                                      | Model of Human Forearm        | This is in a white cardboard box. The label might or might not be facing outwards.                                                                                                                                                                                 |

**Table 1.11.1:** Equipment Needed: Human Forearm

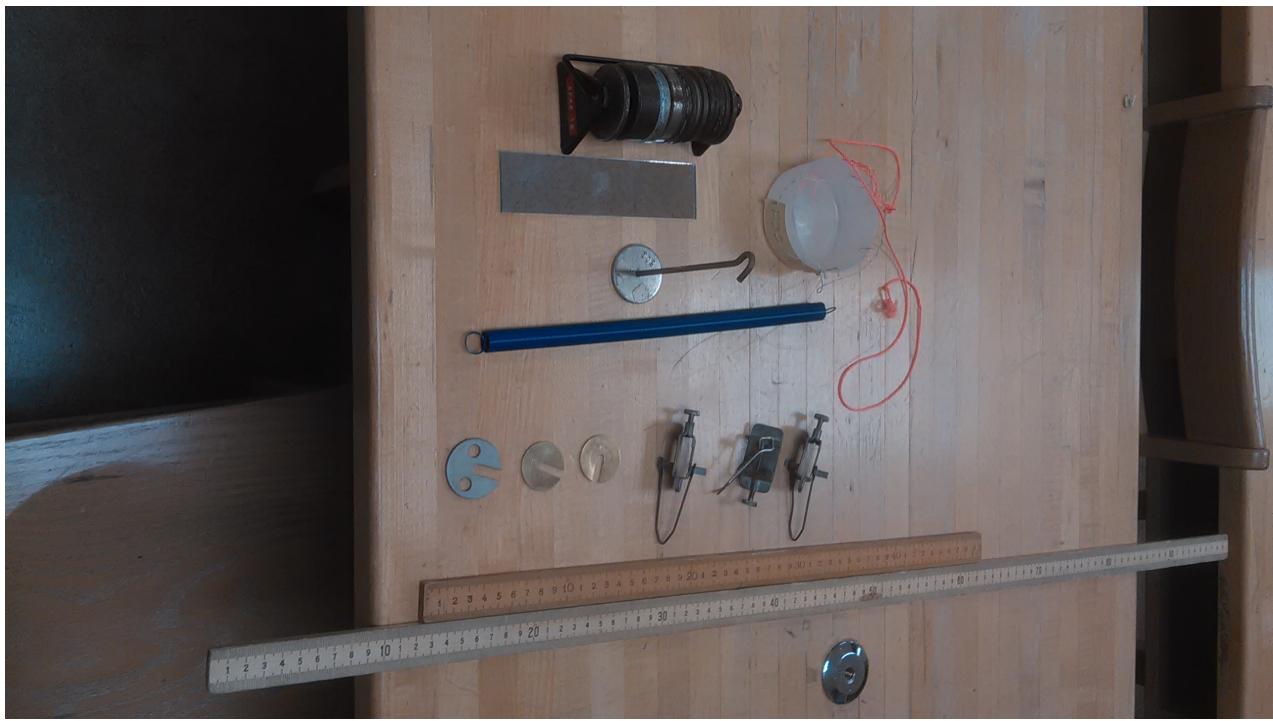


Figure 1.11.2: The necessary equipment laid out.



**Figure 1.11.3:** *If* your instructor wants you to set up the “arm”, then this is how it is connected. Notice that the two clamps on the outer edges hang down and the clamp connected to the spring “hangs” up.



**Figure 1.11.4:** How to set up the support. Recall also [Figure 1.7.3](#) and [Figure 1.7.4](#).

(Updated: November 15, 2017)

A digital version of the lab should be found at <http://physics.thomasmore.edu/Labs/121/c-equilibrium.html>

A PDF version might be found at [equilibrium.pdf](#)

## 1.12 Next Lab

| Location                                  | Equipment | Notes |
|-------------------------------------------|-----------|-------|
| <i>For Each Lab Station</i>               |           |       |
| AF12                                      | 1 track   | .     |
| <i>At the front for students to share</i> |           |       |
| .                                         |           |       |

**Table 1.12.1:** Equipment Needed: [Next Lab](#)

(Updated: October 12, 2017)

A digital version of the lab should be found at <http://physics.thomasmore.edu/Labs/121/c-labname.html>

A PDF version might be found at [labname.pdf](#)



# Class 2

## PHY122L: Elements of Physics (algebra-based, spring)

### 2.1 Skills Review

There is no lab equipment to be set up for this lab. The instructor will provide the students with a worksheet.

(Updated: Jan 8, 2018)

A digital version of the lab should be found at <http://physics.thomasmore.edu/Labs/122/c-122skills.html>

A PDF version might be found at <122skills.pdf>

### 2.2 Archimedes' Principle

You might also review [Archimedes' Principle](#) for another version with less advice for setting up the lab.

| Location                                  | Equipment                                                                                                                                                            | Notes                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     |
|-------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>For Each Lab Station</b>               |                                                                                                                                                                      |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                           |
| AA13                                      | Vernier Caliper                                                                                                                                                      |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                           |
| AC11-12                                   | <i>precise</i> graduated cylinder                                                                                                                                    | You will probably need the tall, skinny ones, but the important thing is to get the cylinders with the finest gradation in their measurement.                                                                                                                                                                                                                                                                                                                                                                                                             |
| AE82                                      | One brass cylinder <b>NOTE:</b> your instructor might want a variety of shapes of objects.                                                                           | These are in a box in the AE82 drawer, but you can use other cylinders from a different drawer or other shapes if indicated by the instructor.                                                                                                                                                                                                                                                                                                                                                                                                            |
| AF35                                      | Balance scale                                                                                                                                                        | The porcelain scales are fragile. Please be careful when handling these balances. These will be placed on the short rods. You should ask your instructor if they want you to do this. There is a picture below.                                                                                                                                                                                                                                                                                                                                           |
| AF44                                      | A tray of masses.                                                                                                                                                    | (I suggest using the cylindrical masses even though they do not need a hanger.) These will be used on the second scale of the scale-balance and should be at least as heavy as the cylinder placed at the table.                                                                                                                                                                                                                                                                                                                                          |
| end of AF                                 | short rod that screws into the table-top.                                                                                                                            | Shorter is better. The balance scale will be placed on top of this and short rods will wobble less.                                                                                                                                                                                                                                                                                                                                                                                                                                                       |
| AG14                                      | Lab Jack                                                                                                                                                             | These need to be tall enough to place the spill container high enough under the balance, which is on the pole, so that the string from the balance reaches into the container. You might need to tie a new string to attach the metal cylinder to the balance.                                                                                                                                                                                                                                                                                            |
| AG44                                      | Spill container                                                                                                                                                      | These are the aluminum "cups" that have a small spout projecting from the side near the top.                                                                                                                                                                                                                                                                                                                                                                                                                                                              |
| <b>At the front for students to share</b> |                                                                                                                                                                      |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                           |
| either AC11 or AG35                       | EITHER the drinking glass (with the slight bulge 3/4 of the way up, in AC11) OR the small beaker (with a pour-spout, in AC11) OR the black and silver cups (in AG35) | This will be used as the "boat" that will carry the pennies. The point is to see how many pennies it takes to sink; the glass makes it difficult to measure volume, but the beaker has a thick lip and a lowered pour-spout. The instructor may want a different boat, such as the aluminum cup in the calorimeter (AG44?) from the <a href="#">specific heat lab</a> ; but that has a very large lip. <b>NOTE:</b> The instructor might want one boat at the front for everybody to use, or the instructor might want each table to have their own boat. |
| either AH18 or AF63                       | a tank for holding several gallons of water                                                                                                                          | There are two cylindrical tanks at AH18; You should only need the smaller tank, but the instructor might want a larger one. There is a fish tank at AF63 that should also work. Your instructor might want you to fill the tank with water. It should be deep enough that the "boat" can sink below the depth of water.                                                                                                                                                                                                                                   |
| AL32                                      | pennies                                                                                                                                                              | You need enough to sink the boat, but each sinking gets the pennies wet. Since the pennies are annoying to dry off, there should be enough pennies to sink at least one (maybe two) boats per group. Before putting the pennies away, you should make sure to dry them off so they don't get moldy. You can use a towel for the most part, but they might have to sit out overnight. <b>NOTE:</b> The instructor might want the pennies distributed among the tables.                                                                                     |


**Table 2.2.1:** Equipment Needed: [Archimedes' Principle](#)



**Figure 2.2.2:** The equipment setup for each group. (Caliper not shown, but needed)



Figure 2.2.3: The equipment setup showing the balance on the pole.

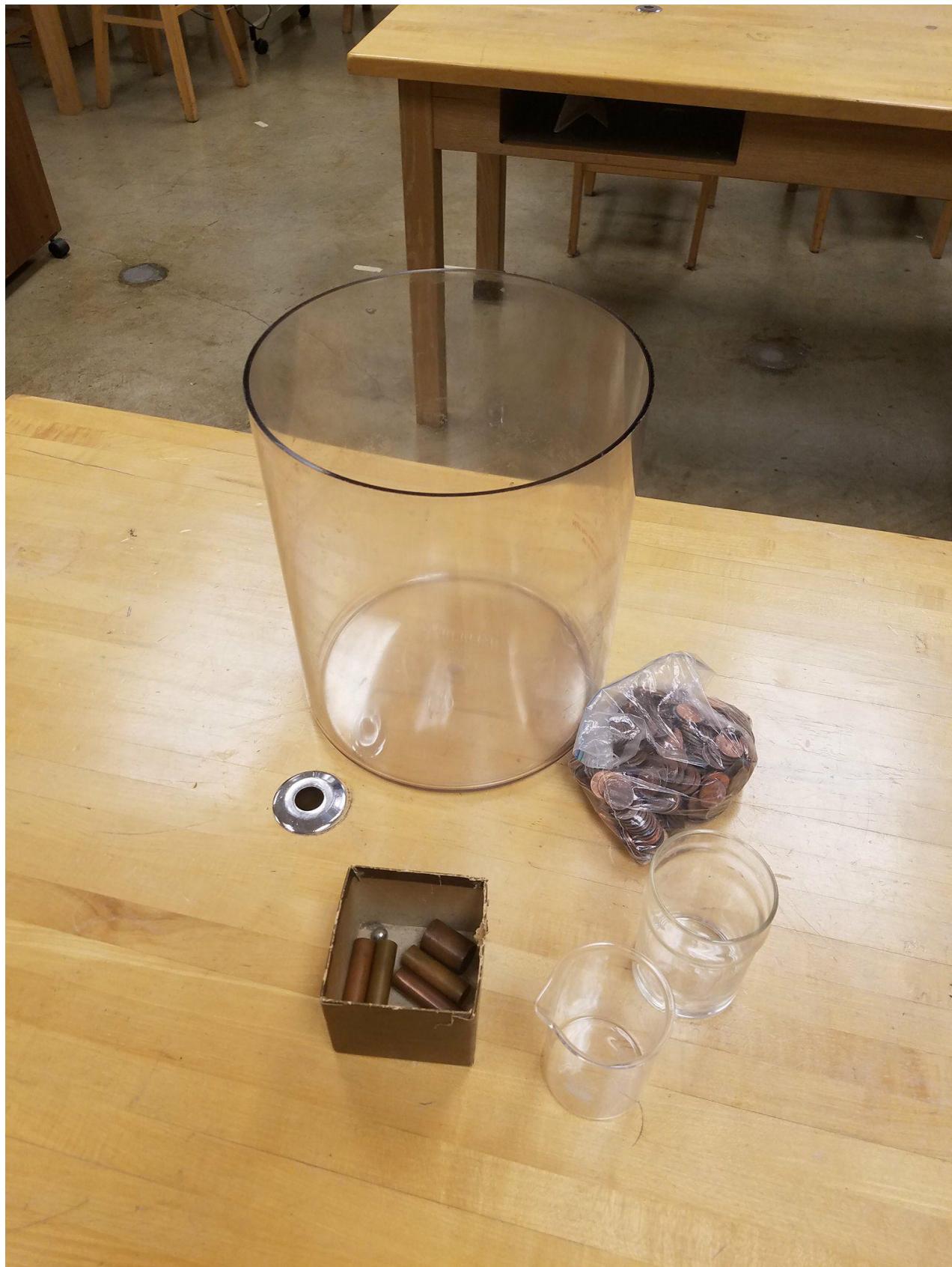

**Figure 2.2.4:** The equipment setup for the front of the room.

(Updated: Jan 9, 2018)

A digital version of the lab should be found at <http://physics.thomasmore.edu/labs/122/c-Archimedes.html>

A PDF version might be found at [Archimedes.pdf](#)

## 2.3 Thermal Expansion

**Alternate Version:** Please use the set-up in [Thermal Expansion](#).

## 2.4 Specific and Latent Heat

**Note:** The faculty member might be doing just specific heat or might be doing both specific and latent heat. Please ask the instructor if the latent heat items are needed.

**Alternate Version:** Your instructor might want the set-up in [Specific Heat Capacity](#).

| Location                                  | Equipment                     | Notes                                                                                                                                                                                                                                                                                                                           |
|-------------------------------------------|-------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>For Each Lab Station</b>               |                               |                                                                                                                                                                                                                                                                                                                                 |
| BB21                                      | thermometer                   | (These are in the drawer under the sink in S224.)                                                                                                                                                                                                                                                                               |
| AC21                                      | 2 beakers                     | One beaker should be 400 mL (for boiling water) and the other should be 1000 mL (for transporting water)                                                                                                                                                                                                                        |
| AD41                                      | 90° clamp                     | This should connect a vertical rod to a horizontal rod                                                                                                                                                                                                                                                                          |
| end of A                                  | 2 small rods                  | One should go into the table hole, the other should cantilever out over the boiling water                                                                                                                                                                                                                                       |
| AE21                                      | string                        | <b>Check with the instructor</b> There are several types of objects that might be used: <ul style="list-style-type: none"> <li>• <i>Most likely</i> AG37 blocks with hooks.</li> <li>• <i>Alternatives:</i> AE82 Spheres with holes, AF44 Brass cylinders, AF44 Aluminum blocks, AF44 some 200-g cylindrical masses.</li> </ul> |
| Various                                   | metal for heating             |                                                                                                                                                                                                                                                                                                                                 |
| AG15                                      | 2 thermal gloves              |                                                                                                                                                                                                                                                                                                                                 |
| AG44                                      | Calorimeter                   | This is the larger aluminum can with a black rubber ring supporting a smaller aluminum can inside. It should have a wooden lid with two holes in the lid. One hole is for the stir-rod (which might be one shelf up) and the other is for the thermometer. One of these holes should have a cork in it.                         |
| AG45                                      | Hot Plate                     |                                                                                                                                                                                                                                                                                                                                 |
| <b>At the front for students to share</b> |                               |                                                                                                                                                                                                                                                                                                                                 |
| AF17?                                     | 2 large cylindrical “buckets” | These are used to hold and insulate ice. The smaller cylinder will have the ice that students use. The larger cylinder will have ice packed around the smaller cylinder. (This way, the outer ice will melt keeping the inner ice cold enough to not melt.)                                                                     |
| Biology Dept.                             | ice                           | They have an ice machine that we can steal ice from. But we need to tell them when we are taking it and we should do it at the beginning of our lab.                                                                                                                                                                            |

**Table 2.4.1:** Equipment Needed: [Specific and Latent Heat](#)



**Figure 2.4.2:** The equipment setup.

(Updated: Feb 2018)

A digital version of the lab should be found at <http://physics.thomasmore.edu/Labs/122/c-spheat.html>

A PDF version might be found at [spheat.pdf](#)

## 2.5 Properties of a Battery

| Location                                  | Equipment                  | Notes                                                                                                                                                                                                                                                                                                                                                                                                                                                                       |
|-------------------------------------------|----------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>For Each Lab Station</b>               |                            |                                                                                                                                                                                                                                                                                                                                                                                                                                                                             |
| AK12                                      | Power Supply               | If this is for the PHY142 class of majors, then use the power supplies at AK55. If this is for the PHY 122 students, then use the power supplies at AK12 and AK22.                                                                                                                                                                                                                                                                                                          |
| AN22                                      | 2 multimeters              | These should all be the same, if possible. Use either the BK2905 or the CM2400T. Sometimes these fuses go bad. In that case, the ammeter portion does not work. It would be best if you could test them before the lab. (You might put 1-2 spares at the front of the room just in case.) When the lab is over, you should ask the instructor if any of the fuses shorted out and replace the fuses in those multimeters.                                                   |
| Electronics Room                          | Resistors                  | <i>This is instructor-dependent.</i> I recommend 8-10 resistors ranging from $300\ \Omega - 800\ \Omega$ (which allows one to stay on the 20 mA setting). I think you can get 2 from each section of the drawers. You should verify (with an ohm-meter) that they are all different resistors. For the PHY 122 students, the instructor might want the resistors laid out in order. For the PHY 142 students, the instructor probably does not want them laid out in order. |
| End Cap of J & K                          | Wires and Wire attachments | Each table should have 3 red and 3 black medium length wires                                                                                                                                                                                                                                                                                                                                                                                                                |
| End Cap of J & K                          | 4 alligator clips          | Check to ensure that the clips form a good connection with the wires and do not slip off easily.                                                                                                                                                                                                                                                                                                                                                                            |
| <b>At the front for students to share</b> |                            |                                                                                                                                                                                                                                                                                                                                                                                                                                                                             |
| AL12 / AL14                               | Batteries                  | There should be a mix of good and bad batteries. Work with the instructor to ensure that if students identify mislabelled batteries, then they should be put into the appropriate bin. The instructor might not want good and bad batteries identified for students, but we still need to keep them separate. There should be a mix of 9 v, 6 v, and 1.5 v. The AA, AAA, C, and D cell batteries are all 1.5 v.                                                             |

**Table 2.5.1:** Equipment Needed: [Properties of a Battery](#)



**Figure 2.5.2:** The equipment set up. This image is for the PHY 142 set-up and does not show all of the equipment described above.

(Updated: Feb 2018)

A digital version of the lab should be found at <http://physics.thomasmore.edu/Labs/122/c-battery.html>

A PDF version might be found at [battery.pdf](#)

## 2.6 Simple Lens

| Location                    | Equipment     | Notes                                                                                                                                                                                                                                                                        |
|-----------------------------|---------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>For Each Lab Station</b> |               |                                                                                                                                                                                                                                                                              |
| AK43                        | Light Source  |                                                                                                                                                                                                                                                                              |
| AK44                        | Optics track  | These are the black tracks that look like the cart-tracks, but they are different. You might note that one side has a ruler attached. This ruler-side should be on the side of the table where the student sits.                                                             |
| AK45                        | Box of lenses | The plastic tub is labeled as having "screens, +10 cm focal length, +20 cm focal length, +25 cm focal length, and -15 cm focal length". Each table should get a screen and one of each of the positive-focal-length lenses. We do not need the negative focal length lenses. |

*At the front for students to share*

**Table 2.6.1:** Equipment Needed: [Simple Lens](#)



**Figure 2.6.2:** The equipment setup.

(Updated: 2018)

A digital version of the lab should be found at <http://physics.thomasmore.edu/Labs/122/c-simplens.html>  
A PDF version might be found at [simplens.pdf](http://physics.thomasmore.edu/Labs/122/c-simplens.pdf)

# Class 3

## NSC 220L: Anything Physics

### 3.1 Meaningful Measurements

Please set this up according to the [Meaningful Measurements](#) from the [PHY 121L: Elements of Physics \(algebra-based, fall\)](#).

### 3.2 Measuring Motion

| Location                                         | Equipment                                | Notes                                                                               |
|--------------------------------------------------|------------------------------------------|-------------------------------------------------------------------------------------|
| <b><i>For Each Lab Station</i></b>               |                                          |                                                                                     |
| AD91                                             | 1 1-meter stick                          | -                                                                                   |
| end of AF                                        | 2-meter stick                            | -                                                                                   |
| AA41 or<br>AA42<br>(drawers)                     | 1 gravity protractor                     | This is the large yellow protractor                                                 |
| <b><i>At the front for students to share</i></b> |                                          |                                                                                     |
| AL32-5                                           | several boxes/bags of pennies            | in sufficient number for each student to have up to 6 pennies                       |
| ?                                                | 2 rolls of masking tape                  | When you find this, tell Dr. Christensen the bin-number so he can update this page. |
| <b><i>set at two tables only</i></b>             |                                          |                                                                                     |
| S224                                             | 1 motion sensor (same as “sonic ranger”) | Should have a black-yellow plug; please connect this to the PASCO interface         |

**Table 3.2.1:** Equipment Needed: [Measuring Motion](#)

(Updated: January 15, 2018)

A digital version of the lab should be found at <http://physics.thomasmore.edu/Labs/220/c-motion.html>

A PDF version of the write-up might be found at [motion.pdf \(128 kB\)](#)

### 3.3 Constant Acceleration

Please set this up according to the [Constant Acceleration](#) from the [PHY 121L: Elements of Physics \(algebra-based, fall\)](#).

### 3.4 Newton's 2<sup>nd</sup> Law on a Linear Track with the Sonic Ranger

Please set this up according to the [Newton's 2<sup>nd</sup> Law on a Linear Track with the Sonic Ranger](#) from the [PHY 121L: Elements of Physics \(algebra-based, fall\)](#).

### 3.5 Dry Sliding Friction

Please set this up according to the [Dry Sliding Friction](#) from the [PHY 121L: Elements of Physics \(algebra-based, fall\)](#).

### 3.6 Hooke's Law

I have not written this yet. Tentatively, please set this up according to the [Hooke's Law / Springs / Peer Review](#) from the [PHY 121L: Elements of Physics \(algebra-based, fall\)](#).

### 3.7 Pendulum / Peer Review

Please set this up according to the [Pendulum / Peer Review](#) from the [PHY 121L: Elements of Physics \(algebra-based, fall\)](#).

### 3.8 Conservation of Energy

Please set this up according to the [Conservation of Energy](#) from the [PHY 121L: Elements of Physics \(algebra-based, fall\)](#).

### 3.9 Bouncing Springs

I have not written this yet. Tentatively, please set this up according to the [Hooke's Law / Springs / Peer Review](#) from the [PHY 121L: Elements of Physics \(algebra-based, fall\)](#).

## Class 4

# PHY141L: General Physics I (calculus-based, spring)

See also [PHY 121L: Elements of Physics \(algebra-based, fall\)](#)



# Class 5

## PHY142L: General Physics II (calculus-based, spring)

See also [PHY122L: Elements of Physics \(algebra-based, spring\)](#)

### 5.1 Ohm's Law (Ohmic material)

*2013 version by Matthew: May need updating..*

| Location                           | Equipment                  | Notes                                 |
|------------------------------------|----------------------------|---------------------------------------|
| <b><i>For Each Lab Station</i></b> |                            |                                       |
| AK55                               | Power Supply               | (JS listed a function generator AK23) |
| End Cap of J & K                   | Wires and Wire attachments |                                       |
| AN22                               | Multimeter                 | Requires 4                            |
| AI25                               | Ammeter                    | Around 3 to 10 amps                   |
| Electronics Room                   | Resisters and Diodes       |                                       |

---

*At the front for students to share*

**Table 5.1.1:** Equipment Needed: [Ohm's Law \(Ohmic material\)](#)



**Figure 5.1.2:** The equipment setup.


(Updated: 2013)

## 5.2 Ohm's Law (non-Ohmic material)

*2013 version by Matthew: May need updating..*

| Location                                  | Equipment                  | Notes      |
|-------------------------------------------|----------------------------|------------|
| <b>For Each Lab Station</b>               |                            |            |
| AK23                                      | Function Generator         |            |
| End Cap of J & K                          | Wires and Wire attachments |            |
| AN22                                      | Multimeter                 | Requires 4 |
| Electronics Room                          | Resisters and Diodes       |            |
| <b>At the front for students to share</b> |                            |            |

**Table 5.2.1:** Equipment Needed: [Ohm's Law \(non-Ohmic material\)](#)

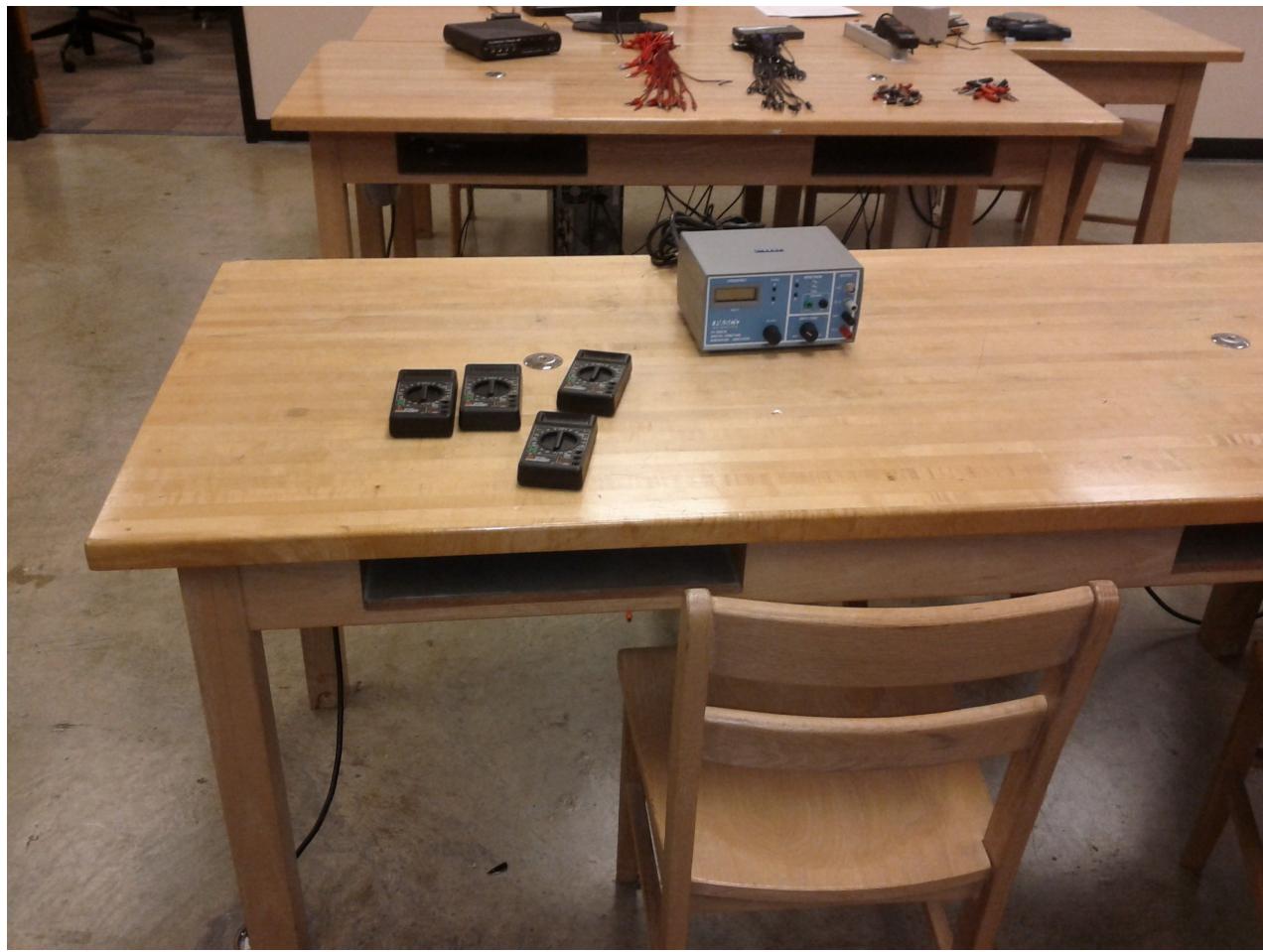

**Figure 5.2.2:** The equipment setup.

(Updated: 2013)

### 5.3 Mapping the Electric Field

*2013 version by Matthew: May need updating..*

| Location                    | Equipment                              | Notes                                                                                                                                                                                                        |
|-----------------------------|----------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>For Each Lab Station</b> |                                        |                                                                                                                                                                                                              |
| AK55                        | Power Supply                           |                                                                                                                                                                                                              |
| End cap of J & K            | Wire and Wire Supplies                 | JS also lists                                                                                                                                                                                                |
| AN22                        | Multimeter                             | <ul style="list-style-type: none"> <li>• AJ43</li> <li>• AJ25: Small Black ones &amp; Large Blue ones</li> <li>• AJ24: Small Blue ones &amp; Boxed ones;</li> <li>• AI15: Ones with three buttons</li> </ul> |
| AG24                        | Field Plate Apparatus and Field Plates |                                                                                                                                                                                                              |

*At the front for students to share*

**Table 5.3.1:** Equipment Needed: [Mapping the Electric Field](#)



**Figure 5.3.2:** The equipment setup.

(Updated: 2013)

## 5.4 Resistor Bridge

*2013 version by Matthew: May need updating..*

| Location                                  | Equipment         | Notes |
|-------------------------------------------|-------------------|-------|
| <b>For Each Lab Station</b>               |                   |       |
| AJ43                                      | Bridge Apparatus  |       |
| AJ24                                      | Resistor Coil Set |       |
| AJ23                                      | Resistor Box      |       |
| AK23                                      | Power Supply      |       |
| End of Isle J & K                         | Wires & Supplies  |       |
| AI15                                      | Galvanometer      |       |
| <i>At the front for students to share</i> |                   |       |

**Table 5.4.1:** Equipment Needed: Resistor Bridge



**Figure 5.4.2:** The equipment set up.


(Updated: 2013)

## 5.5 Magnetic Field Strength

*2013 version by Matthew: May need updating..*

| Location                                         | Equipment        | Notes |
|--------------------------------------------------|------------------|-------|
| <b><i>For Each Lab Station</i></b>               |                  |       |
| AM36                                             | Magnet           |       |
| AG14                                             | 2 Lab Jacks      |       |
| AF34                                             | Ruler            |       |
| AF15                                             | 2 Blocks of Wood |       |
| Room S224 Closet                                 | Magnetic Sensor  |       |
| <b><i>At the front for students to share</i></b> |                  |       |

**Table 5.5.1:** Equipment Needed: [Magnetic Field Strength](#)

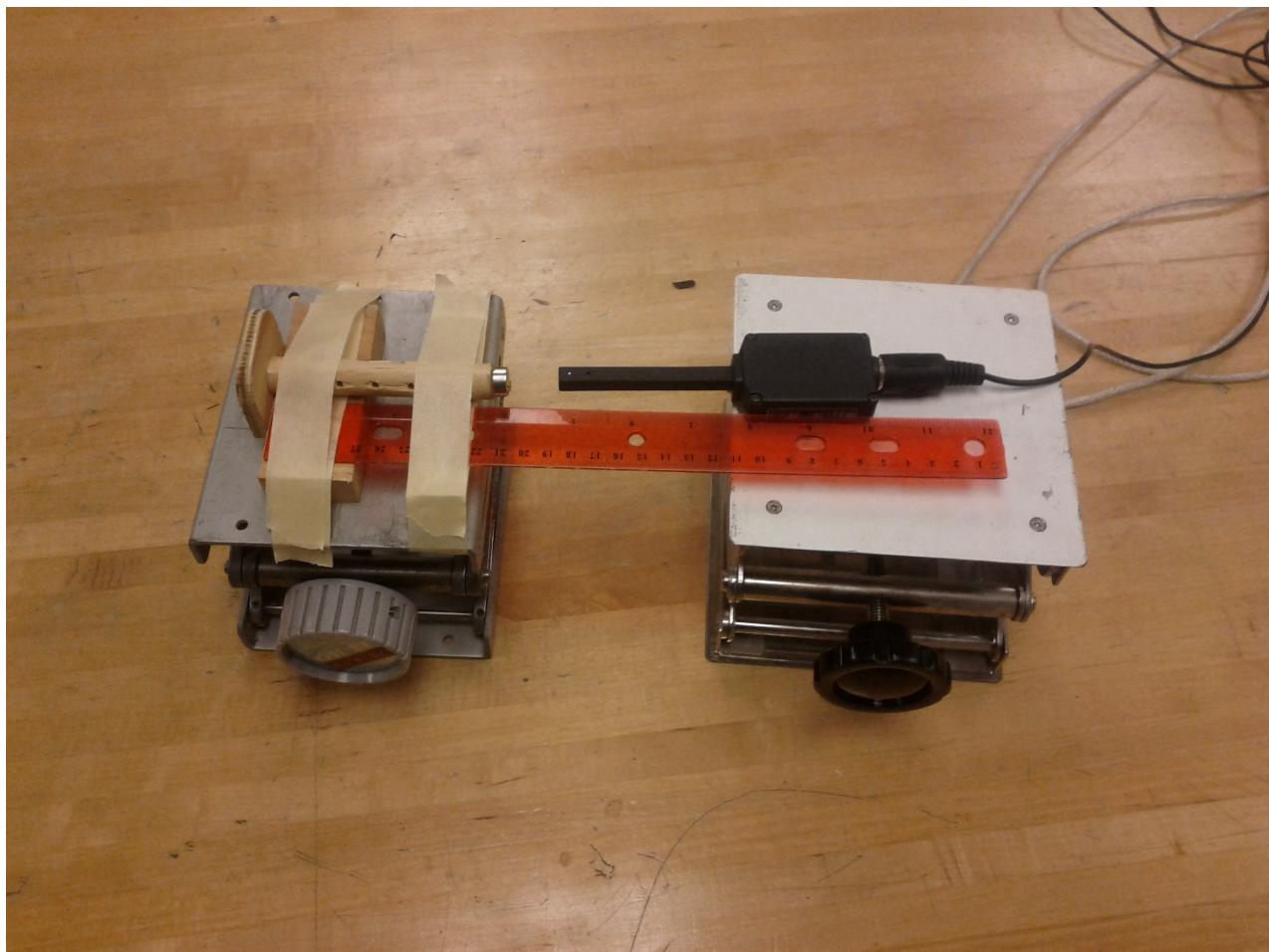

**Figure 5.5.2:** The equipment setup.

(Updated: 2013)

## 5.6 High Current Force

*2013 version by Matthew: May need updating..*

| Location                    | Equipment                 | Notes                  |
|-----------------------------|---------------------------|------------------------|
| <i>For Each Lab Station</i> |                           |                        |
| AG14                        | Lab Jack                  |                        |
| AH12                        | Laser                     |                        |
| BA14                        | Current apparatus         |                        |
| AE61                        | Double-throw Switch       |                        |
| AF34                        | Ruler                     |                        |
| End of Isle J &K            | Wire and Supplies         |                        |
| AK12                        | High Current Power Supply |                        |
| AN22                        | Multimeter                |                        |
| AA13                        | Caliper                   | (inconsistent listing) |
| AF43                        | Small mass set            | (inconsistent listing) |

---

*At the front for students to share*

---

**Table 5.6.1:** Equipment Needed: [High Current Force](#)**Figure 5.6.2:** The equipment setup.

(Updated: 2013)

## 5.7 Magnet in a Tube

*2013 version by Matthew: May need updating..*

| Location                                         | Equipment               | Notes |
|--------------------------------------------------|-------------------------|-------|
| <b><i>For Each Lab Station</i></b>               |                         |       |
| AH63                                             | Magnetic pipe apparatus |       |
| white box                                        | Magnet and non-magnet   |       |
| -                                                | Foam Pads               |       |
| S224 Closet                                      | Force Sensor            |       |
| <b><i>At the front for students to share</i></b> |                         |       |

**Table 5.7.1:** Equipment Needed: [Magnet in a Tube](#)

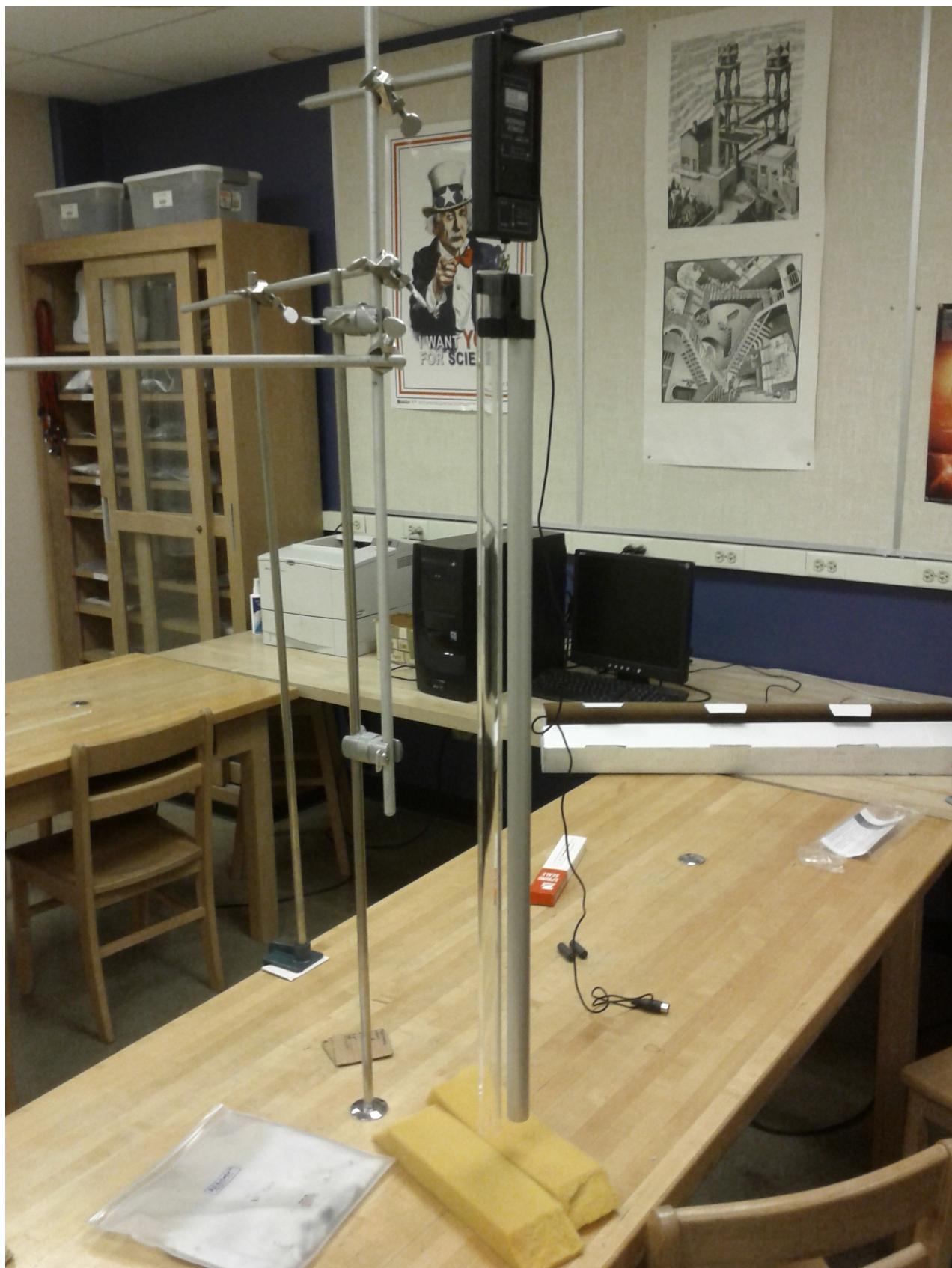


Figure 5.7.2: The equipment set up.


(Updated: 2013)

## 5.8 Capacitor and Variable Resistor

*2013 version by Matthew: May need updating..*

| Location                                  | Equipment              | Notes |
|-------------------------------------------|------------------------|-------|
| <b>For Each Lab Station</b>               |                        |       |
| AK23                                      | Function Generator     |       |
| End of Isle J & K                         | Wire and Supplies      |       |
| AN22                                      | Multimeter             |       |
| AJ13                                      | Resistor Apparatus     |       |
| AJ26                                      | Varying Resistor Board |       |
| <i>At the front for students to share</i> |                        |       |

**Table 5.8.1:** Equipment Needed: Capacitor and Variable Resistor

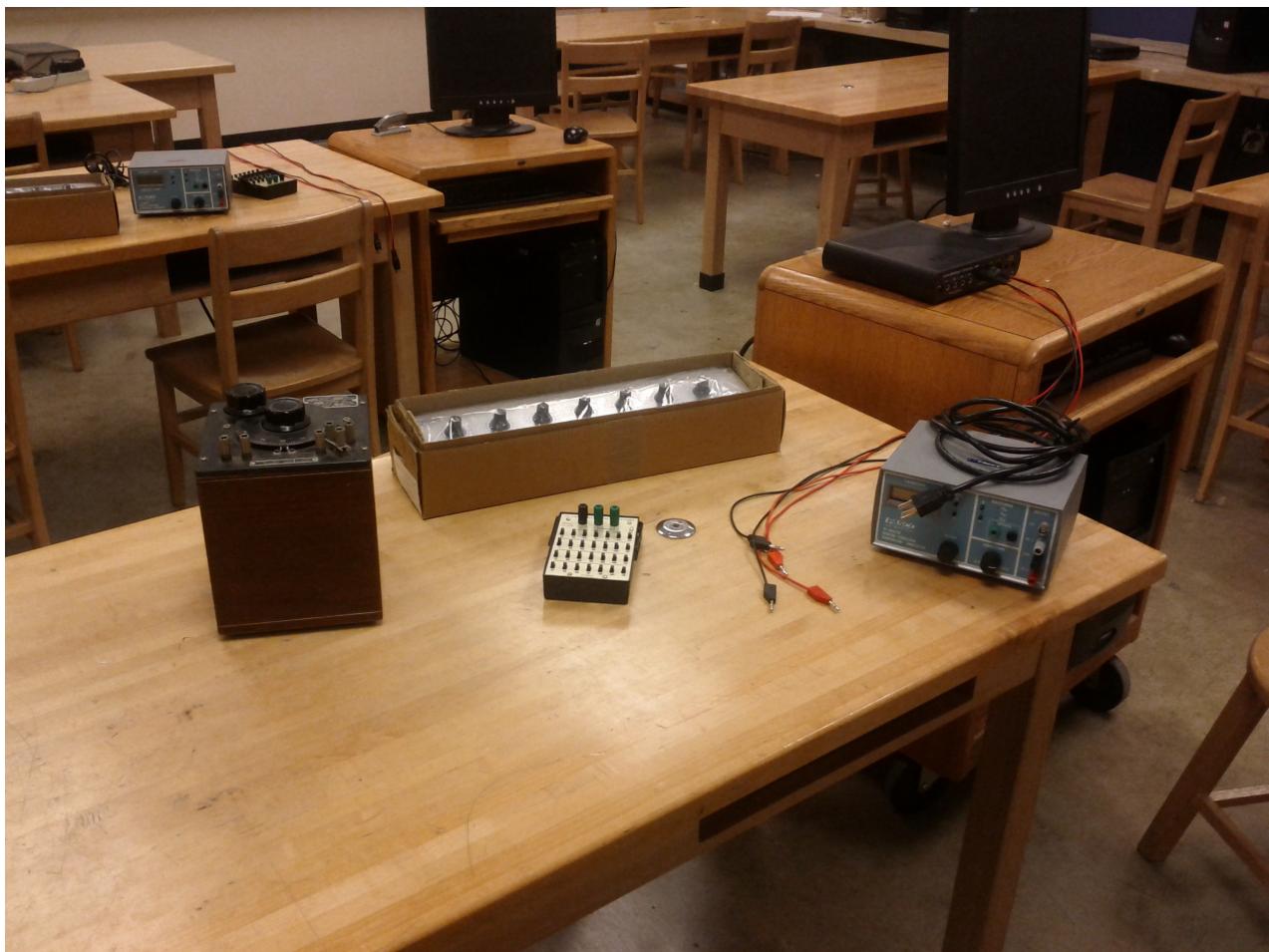

**Figure 5.8.2:** The equipment setup.

(Updated: 2013)

# Class 6

## PHY241L: General Physics III (calculus-based, spring)

See also [PHY122L: Elements of Physics \(algebra-based, spring\)](#)

### 6.1 Archimedes' Principle

*2013 version by Matthew: May need updating..* Please also review [Archimedes' Principle](#) for another version with additional advice for setting up the lab.

| Location                                  | Equipment             | Notes   |
|-------------------------------------------|-----------------------|---------|
| <i>For Each Lab Station</i>               |                       |         |
| AA13                                      | Vernier Caliper       |         |
| AC11 (aisle E)                            | Beaker                |         |
| AE82                                      | Objects               |         |
| AG44                                      | Overflow Cup          |         |
| AF35                                      | Scales                |         |
| AG14                                      | Lab Jack              |         |
| AG35                                      | Black and silver cups | "Boats" |
| AL32                                      | Pennies               |         |
| <i>At the front for students to share</i> |                       |         |
| AF63                                      | Fish Tank             |         |

**Table 6.1.1:** Equipment Needed: [Archimedes' Principle](#)



Figure 6.1.2: The equipment setup.



Figure 6.1.3: The equipment set up.

(Updated: 2013)

## 6.2 Standing Waves

*2013 version by Matthew: May need updating..*

| Location                                  | Equipment          | Notes |
|-------------------------------------------|--------------------|-------|
| <i>For Each Lab Station</i>               |                    |       |
| AK23                                      | Function generator |       |
| End of isle F                             | two-meter stick    |       |
| AG14                                      | Lab jack           |       |
| AG16                                      | Lead Brick         |       |
| AD72 & AE81                               | Pulley system      |       |
| AF13                                      | Vibrator           |       |
| Inside of the vibrator box                | String             |       |
| AF44                                      | Masses and Hooks   |       |
| End cap of isle J & K                     | Wires              |       |
| <i>At the front for students to share</i> |                    |       |
| AN24                                      | Strobe light       |       |

**Table 6.2.1:** Equipment Needed: [Standing Waves](#)



**Figure 6.2.2:** The equipment setup.


(Updated: 2013)

### 6.3 Speed of Sound in a Tube

*2013 version by Matthew: May need updating..*

| Location                                  | Equipment          | Notes                                                    |
|-------------------------------------------|--------------------|----------------------------------------------------------|
| <b>For Each Lab Station</b>               |                    |                                                          |
| AK63                                      | Speaker            |                                                          |
| Room S224 Closet                          | Sound Sensor       |                                                          |
| AF12                                      | Cardboard Tube     | Use the cardboard tubes used to hold the aluminum tracks |
| End of isle F                             | two-meter stick    |                                                          |
| AK23                                      | Function Generator |                                                          |
| <b>At the front for students to share</b> |                    |                                                          |

**Table 6.3.1:** Equipment Needed: [Speed of Sound in a Tube](#)

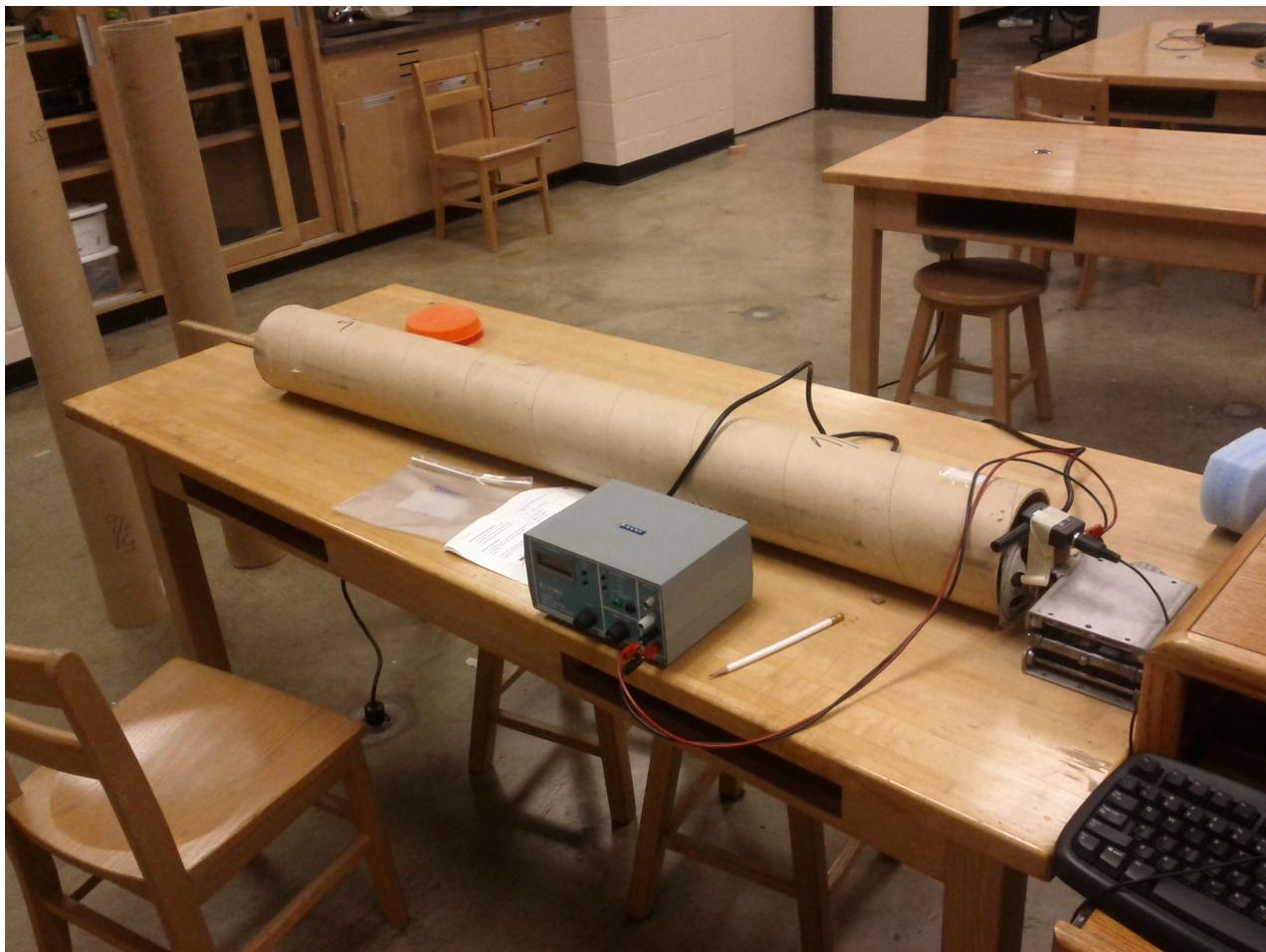

**Figure 6.3.2:** The equipment setup.

## 6.4 Thermal Expansion

| Location                                  | Equipment                                | Notes                                                                                                                                                                                                                                 |
|-------------------------------------------|------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>For Each Lab Station</b>               |                                          |                                                                                                                                                                                                                                       |
| AE41                                      | Stoppers                                 |                                                                                                                                                                                                                                       |
| AF (end)                                  | Two meter-sticks                         | One at the table to measure the rods and one taped to the wall with masking tape                                                                                                                                                      |
| AF31?                                     | Glass jar with aluminum cap with a spout |                                                                                                                                                                                                                                       |
| AF34                                      | Ruler                                    |                                                                                                                                                                                                                                       |
| AG14                                      | Lab jack                                 | Set the lab jack on top of a short lab stool that rests on the lab table. The laser will be placed on top of this during the lab. (The laser should be on the table until the lab starts so it is less likely to fall to the ground.) |
| AG15                                      | Thermal Gloves                           |                                                                                                                                                                                                                                       |
| AG45                                      | Hot plate                                | We might only have 5                                                                                                                                                                                                                  |
| AG45                                      | Tray                                     | to catch the spilled water                                                                                                                                                                                                            |
| AG54                                      | Box of rubber tubing                     | You should only use the plastic tubes, because the rubber ones can melt                                                                                                                                                               |
| AG55                                      | Base for linear expansion                |                                                                                                                                                                                                                                       |
| AG55                                      | Box of black corks                       | These probably only fit the central hole (for the thermometer)                                                                                                                                                                        |
| AG55                                      | Towel                                    |                                                                                                                                                                                                                                       |
| AG56                                      | Metal tube                               | This has small spouts on either end and an opening in the middle. <b>Be sure that the tube has a cork in each end</b>                                                                                                                 |
| AH12                                      | Laser                                    |                                                                                                                                                                                                                                       |
| BA11 (in S224)                            | Tripod Mirror                            |                                                                                                                                                                                                                                       |
| BB21 (in S224)                            | Thermometer                              |                                                                                                                                                                                                                                       |
| <b>At the front for students to share</b> |                                          |                                                                                                                                                                                                                                       |
| AG56                                      | Box of expansion rods                    | There should be aluminum, brass, and steel rods. They should all be the same length and none of them should be hollow.                                                                                                                |
| Pegboard                                  | Roll of masking tape                     |                                                                                                                                                                                                                                       |
| Toolbox                                   | Needle nose pliers                       | In case a hot rod needs to be adjusted                                                                                                                                                                                                |

**Table 6.4.1:** Equipment Needed: [Thermal Expansion](#)



**Figure 6.4.2:** The equipment setup.

(Updated: Feb, 2018)

## 6.5 Boyle's Law

*2013 version by Matthew: May need updating..*

| Location                                  | Equipment             | Notes                                          |
|-------------------------------------------|-----------------------|------------------------------------------------|
| <b>For Each Lab Station</b>               |                       |                                                |
| End of isle G & H                         | Boyle's law apparatus |                                                |
| Faucet                                    | water                 | (fill the apparatus to the appropriate height) |
| <b>At the front for students to share</b> |                       |                                                |

**Table 6.5.1:** Equipment Needed: Boyle's Law



Figure 6.5.2: The equipment setup.

(Updated: 2013)

## 6.6 Specific Heat Capacity

*2013 version by Matthew: May need updating..*

**Alternate Version:** Your instructor might want the set-up in [Specific and Latent Heat](#).

| Location                                         | Equipment      | Notes             |
|--------------------------------------------------|----------------|-------------------|
| <b><i>For Each Lab Station</i></b>               |                |                   |
| AC11                                             | Beaker         |                   |
| AG44                                             | Calorimeter    | the insulated cup |
| AG45                                             | Hot Plate      | only have 5       |
| Room S224, BB21                                  | 2 Thermometers |                   |
| End of isle F                                    | Rods           |                   |
| AG37                                             | Objects        |                   |
| AF15 (AG15?)                                     | Gloves         |                   |
| AE21                                             | string         |                   |
| <b><i>At the front for students to share</i></b> |                |                   |

**Table 6.6.1:** Equipment Needed: [Specific Heat Capacity](#)



**Figure 6.6.2:** The equipment set up (revised).



**Figure 6.6.3:** The equipment setup (old version, but better view of most equipment).

(Updated: 2013)

## 6.7 Newton's Law of Cooling

*2013 version by Matthew: May need updating..*

| Location                                  | Equipment                       | Notes |
|-------------------------------------------|---------------------------------|-------|
| <i>For Each Lab Station</i>               |                                 |       |
| AG45                                      | Hot Plate                       |       |
| AG11                                      | Large and Small Beaker          |       |
| End of isle F                             | Stand to hold objects and Wires |       |
| Room S224, BB21                           | Thermometer                     |       |
| AF15                                      | Gloves                          |       |
| Either in lab room or room S224           | Computer with Pasco             |       |
| AG58                                      | Bottles                         |       |
| <i>At the front for students to share</i> |                                 |       |

**Table 6.7.1:** Equipment Needed: Newton's Law of Cooling



**Figure 6.7.2:** The equipment setup.

(Updated: 2013)


## 6.8 Snell's Law of Refraction

*Different instructors have different preferences.* some will want the lasers, others will not.

| Location                    | Equipment                             | Notes                                                                                                                                        |
|-----------------------------|---------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------|
| <b>For Each Lab Station</b> |                                       |                                                                                                                                              |
| AH34                        | Thin Refraction squares and triangles |                                                                                                                                              |
| AH44                        | Cardboard Plates                      |                                                                                                                                              |
| AA43                        | Pins                                  |                                                                                                                                              |
| AA42                        | Protractor                            |                                                                                                                                              |
| AH44                        | Red plates                            | <b>Option:</b> If the instructor wants <i>reflection and refraction</i> , then they <i>might</i> want to use this as the reflecting surface. |
| AF43                        | plane mirror                          | <b>Option:</b> If the instructor wants <i>reflection and refraction</i> , then they <i>might</i> want to use this as the reflecting surface. |
| AH12                        | Laser                                 | <b>Note:</b> Not all instructors want the laser; please ask.                                                                                 |
| AG14                        | Lab jack                              | <b>Note:</b> Only use this if using the laser.                                                                                               |

*At the front for students to share*

**Table 6.8.1:** Equipment Needed: [Snell's Law of Refraction](#)

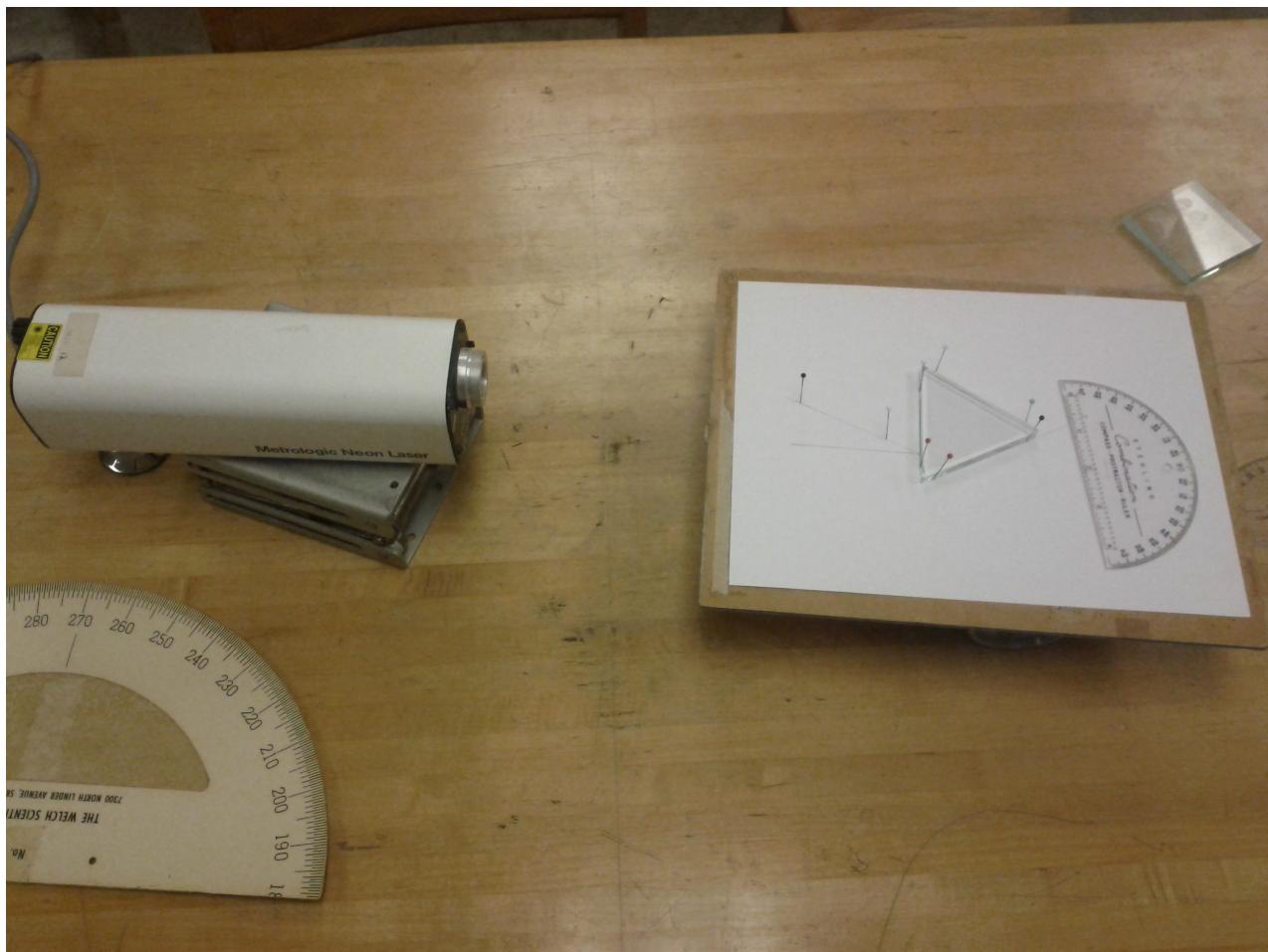

**Figure 6.8.2:** The equipment setup.



**Figure 6.8.3:** The equipment setup; version 2.



**Figure 6.8.4:** A close-up view of why you need two jacks (to get the laser to line up with the side of the lens). If the instructor does not want the lasers, then you also do not need the lab jacks.

(Updated: 2013)

## 6.9 Reflection and Refraction

*2013 version by Matthew: May need updating..*

| Location                    | Equipment         | Notes |
|-----------------------------|-------------------|-------|
| <b>For Each Lab Station</b> |                   |       |
| AH34                        | Refraction Plates |       |
| AH45                        | Cardboard Plates  |       |
| AH45                        | Red Plastic Lens  |       |
| AH12                        | Lamps             |       |
| AG14                        | Lab Jack          |       |
| AA21                        | Hinged Mirror     |       |
| AA43                        | Pins              |       |
| AA42                        | Protractor        |       |
| AG12                        | Light Pipe        |       |

*At the front for students to share*

**Table 6.9.1:** Equipment Needed: [Reflection and Refraction](#)



**Figure 6.9.2:** The equipment setup.

## 6.10 Radiation

*2013 version by Matthew: May need updating..*

| Location                                         | Equipment           | Notes |
|--------------------------------------------------|---------------------|-------|
| <b><i>For Each Lab Station</i></b>               |                     |       |
| (ASK)                                            | Radioactive Samples |       |
| AG16                                             | Lead Bricks         |       |
| -                                                | Stop Watch          |       |
| AM53                                             | Geiger counter      |       |
| AM43                                             | Shields             |       |
| <b><i>At the front for students to share</i></b> |                     |       |

Table 6.10.1: Equipment Needed: [Radiation](#)

(Updated: 2013)

## 6.11 Interferometer

*2013 version by Matthew: May need updating..*

| Location                                         | Equipment      | Notes |
|--------------------------------------------------|----------------|-------|
| <b><i>For Each Lab Station</i></b>               |                |       |
| AN32                                             | Interferometer |       |
| AH12                                             | Laser          |       |
| AK4                                              | Screen         |       |
| AG14                                             | Lab Jacks      |       |
| AG12                                             | Hand Pump      |       |
| <b><i>At the front for students to share</i></b> |                |       |

Table 6.11.1: Equipment Needed: [Interferometer](#)



**Figure 6.11.2:** The equipment setup.

(Updated: 2013)